# MEDISANA®

## TargetScale® 2

Körperanalysewaage  
Body Analysis Scale

- DE** Kurzanleitung
- GB** Quick Start Guide
- FR** Notice abregée
- IT** Guida di avvio rapido
- ES** Guía Rápida
- NL** Beknopte handleiding
- RU** Быстрый пуск

CE

Art. No. 40419

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## IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- Use the device only according to its intended purpose as specified in the instruction manual.
- The warranty will be invalidated if the device is used for purposes other than those for which it is intended.
- This device is not suitable for pregnant women!
- Weighing only is no problem, but bodyfat analysis should be avoided for children because of the current impulse. This impulse is flowing through the body during the bodyfat measurement process.
- These scales are not suitable for people with pacemakers or other medical implants. The body fat information may be inaccurate in people with diabetes or other medical/physical limitations. The same also applies to very highly trained athletes.
- Do not use the device if it is not working correctly, if it has been dropped or has fallen into water or has been damaged.

Protect the unit against moisture. Keep the scale away from water. The surface of the scales is slippery when wet. Keep the surface dry.

### Risk of slipping!


Do not stand on the scale with wet feet. Do not stand on the scale wearing socks.


### Risk of tipping over!

Do not place the scale on an uneven surface. When weighing yourself, do not stand on one side or the edge of the scales.

## Simple body weight measurement using the "Step-On" feature

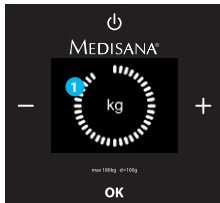
1 During the measurement remain still until the circle is complete.

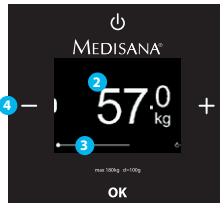
- 2 Your weight is being displayed.
- 3 The timeline displays the result for 5 seconds.
- 4 Using the minus touch key the timeline can be reset. The result will be displayed again.

## Set up personal details on the scale

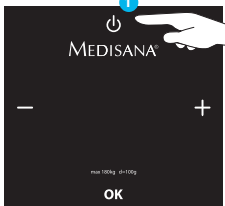1 Tap the ON/OFF touch key briefly.

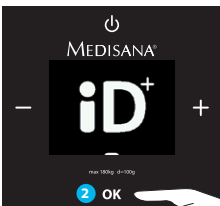2 Tap the OK touch key.


- 3 Select a new avatar (user icon) using the +/- touch keys.
- 4 Confirm your selection using the OK touch key. Follow the additional steps on the display.

### Pair your TargetScale with your iPhone®/iPad®/iPod touch®


The scale and the iPhone®/iPad®/iPod touch® have to be paired once using Bluetooth. Please install the latest VitaDock app and perform the following steps.

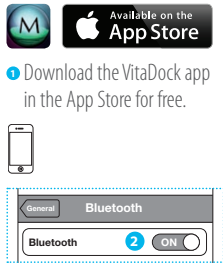- 1 Download the VitaDock app in the App Store for free.
- 2 Activate Bluetooth on your iOS device.


- 3 Tap and hold the ON/OFF touch key until...
- 4 ...the Bluetooth symbol flashes on the display.

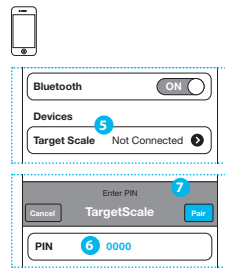- 5 Select TargetScale.
- 6 Enter 0000 as PIN.
- 7 Tap Pair.


- 8 The VitaDock wizard starts after 5 seconds. Allow the devices to communicate.

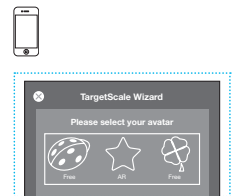
- 9 Please perform the settings in the VitaDock installation assistant (wizard) to benefit from all features of your TargetScale. In case the wizard does not start automatically, you can find it in the TargetScale settings in the app.

- 10 You can now start the measurement.

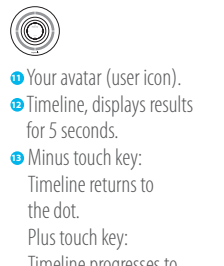
- 11 Your avatar (user icon).
- 12 Timeline, displays results for 5 seconds.
- 13 Minus touch key: Timeline returns to the dot.
- 14 Plus touch key: Timeline progresses to the next dot.

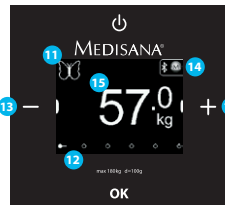
- 14 Bluetooth icon constantly on: Connection between scale and iPhone.
- VitaDock icon constantly on: Connection between scale and VitaDock app.

- 15 The following results will be shown on the display:
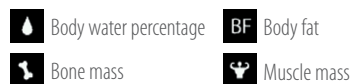
### Delete and reset user settings

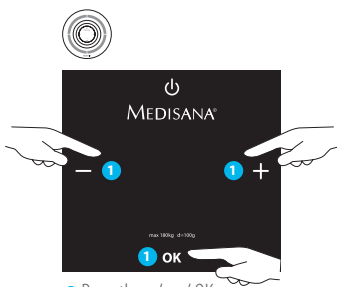- 1 Press the - / + / OK sensor buttons at the same time.

- 2 Tap and hold the touch keys until the circle around the dustbin icon is complete.

- 3 The ID icon is being displayed above the dustbin icon. Now you have the following options:
  - Option 1: Delete single user icon: Tap the OK touch key.
  - Option 2: Reset to factory settings: Tap the + touch key.

- Option 1: Delete single user symbols: All active user icon are being displayed. Select the user icon to be deleted by using the +/- touch key. Tap and hold the OK touch key until the circle around the user icon is complete. All user settings and measurement data of this user will be deleted.

- Option 2: Reset to factory settings: The factory icon is being displayed. Tap and hold the OK touch key until the circle around the factory icon is complete.

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