


QR-7024

Welcome to BarberBoss and congratulations on your new purchase! Secure your peace of mind by activating your 12-month warranty and claiming an extra 6 months for free at www.barberboss.co.uk/barberboss-product-warranty. Beyond our warranty, you are also covered by your country's legal consumer rights, ensuring complete protection for your purchase.

Product Features

- **Advanced 360° Rotating Tip with Gravity Ball Design:** Floss easily in any direction, even upside down, for thorough cleaning. With low-frequency, high-water volume technology, the BarberBoss water flosser effectively removes food particles and plaque while protecting your gums and promoting oral health.
- **Generous 300ml Water Tank:** Large capacity tank allows for uninterrupted use, providing enough water for a full session without needing to refill.
- **5 cleaning mode with long-lasting battery:** With five adjustable cleaning modes, you can customise your oral care to suit your needs for a personalised clean every time. The powerful 2500mAh battery ensures long-lasting performance, providing 20-30 days of use on a single charge for hassle-free convenience.
- **IPX7 Waterproof:** The water flosser is designed with an IPX7 waterproof rating, ensuring all key parts like the motor, battery, and circuit board are safely protected from water damage. This makes it durable and safe to use, even in wet environments.

Product Structure and Accessories

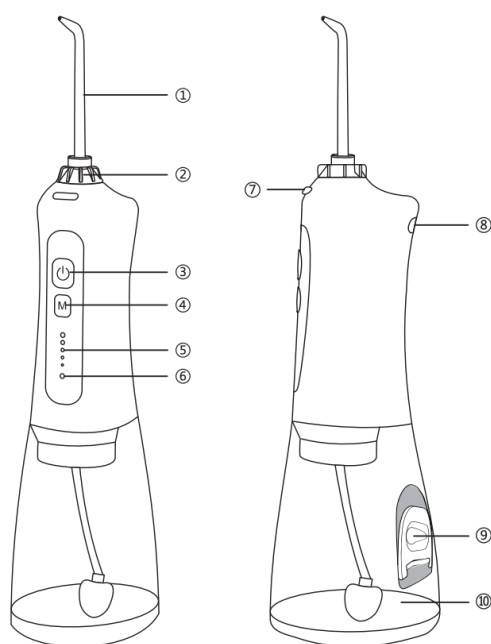

1. Nozzle tip
2. Tip rotation dial
3. Power Button
4. Mode selection button

5. mode indicator lights
6. charging indicator light
7. Tip release button
8. Charging port
9. Water tank cap
10. Water tank




1



2



3

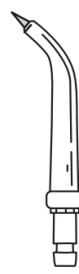

4



5



6

1. USB Charging Cable
2. Tongue Scraper Tip
3. Orthodontic Tip
4. Periodontal Tip
5. Standard Tip
6. Plaque Seeker Tip

Product Details

Product Name	QR-7024 Oral Irrigator
Model Number	QR-7024
Charging Input	5VDC 1.0A
Water Tank Capacity	300ml
Charging Time	5 Hours
Usage Time	20-30 Days

What's Included

Water Flosser	x1
Standard Tip	x2
Orthodontic Tip	x1
Periodontal Tip	x1

Tongue Scraper Tip	x1
Plaque Seeker Tip	x1
Charging Cable	x1

Nozzle Tips: Function and Usage

Nozzle type	Usage	Recommended Replacement Interval
Standard Tip	Ideal for everyday cleaning and removing food debris between teeth and gums.	Every 6 months
Orthodontic Tip	Designed for cleaning around braces, brackets, and other orthodontic devices.	Every 3 months
Periodontal Tip	Gently cleans deep below the gumline to treat gum disease and improve health.	Every 3 months
Tongue Cleaner Tip	Helps remove bacteria and debris from the tongue, promoting fresher breath.	Every 6 months
Plaque Seeker Tip	Targets and removes hard-to-reach plaque around implants, crowns, and bridges.	Every 6 months

Instructions For Use

- Power On/Off:
 - Start Flossing:** Insert the tip, then press the power button to begin irrigating. The device has a 2-minute auto-stop timer or can be turned off manually by pressing the power button again.
 - Turn Off:** Once finished, press the power button to turn off the device.
 - Empty and Drain:** Empty any remaining water from the tank, then run the device for 2-3 seconds to drain residual water.
- Inserting and removing the Nozzle
 - Insert the nozzle into the top of the water flosser handle, pressing firmly until you hear a click, confirming it is securely in place (Figure 1).
 - To remove the nozzle, press the tip release button and gently pull the nozzle out.
 - Note:**
 - Do not press the tip release button while the device is operating.
 - If the nozzle head doesn't spray water immediately after turning on the device (particularly after first use or long periods of inactivity), follow these

steps: 1. Add a small amount of water to the tip insertion area (Figure 2). 2. Insert a new tip and press the power button to activate the device.

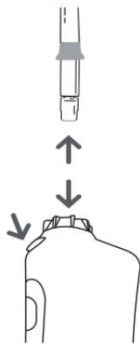


Figure 1

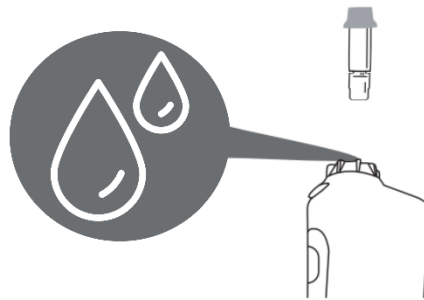


Figure 2

3. Installing and Filling the Water Tank:

- Open the water tank cap, fill with water, and close the cap securely (Figure 3).
- To detach or attach the water tank, simply twist it.
 - Place the device on a flat surface.
 - Hold the water tank with your left hand and twist the upper part to the right with your right hand to detach.
 - To reattach, twist in the opposite direction (Figure 4).
- Note: Use water or a dental professional-recommended solution only when filling the tank.

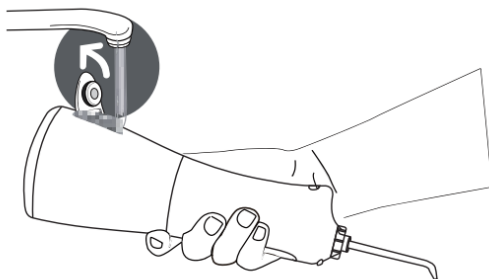


Figure 3

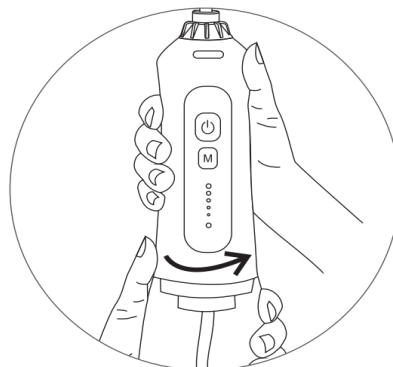


Figure 4

4. Pressure Mode settings

This water flosser offers 5 pressure modes, easily adjustable by pressing the “M” button. The pressure settings are as follows:

- **Highest (110 PSI):** Powerful clean
- **Strong:** Deep cleaning
- **Normal:** Daily cleaning
- **Gentle:** For sensitive teeth
- **Low (40 PSI):** Ideal for beginners

5. How To Charge

- Connect the provided USB cable directly to the device. Then plug the USB into a 5V1A USB adapter. The product must be charged via a laptop/computer/power bank USB port (Figure 5). After charging remove the power plug from the USB outlet and then remove the USB out of the product.
- Charge the product fully before you use it for the first time and after a long period of disuse.
- Make sure the device is switched OFF before charging.
- The LED light will light up, confirming the charging process has begun.
- A full charge takes 3-5 hours, providing up to 30 days of operation. For optimal performance, it's best to use the device when fully charged.

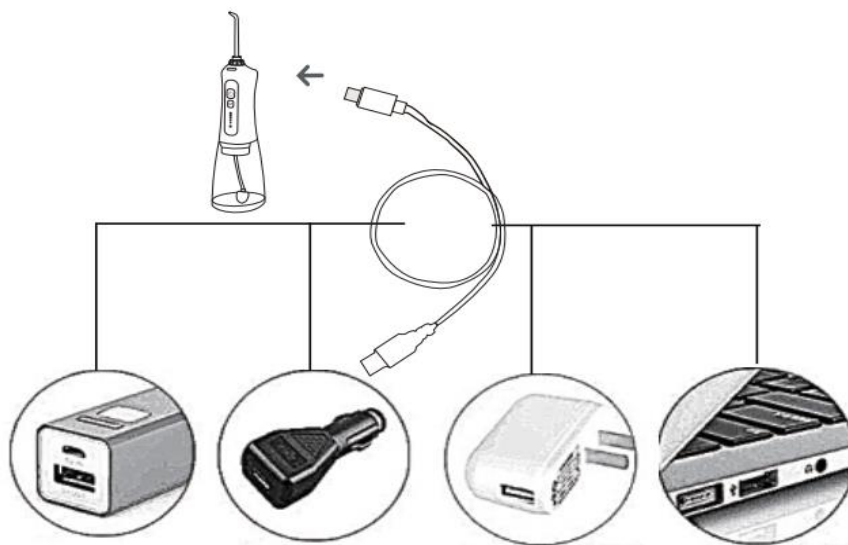


Figure 5

Cleaning & Maintenance

Please note: Although this appliance is waterproof (IPX7 Rating), do not immerse the device in water or any other liquid. To maintain optimal performance and longevity, make sure that the product is cleaned and thoroughly dried before its next use.

- Before cleaning, ensure the device is switched OFF and disconnected from any charging cables.
- Empty the water tank after each use and run the device for 2-3 seconds to drain any remaining water to prevent bacteria buildup.
- Always use only water or neutral cleaning agents to clean the appliance. Never use caustic or abrasive agents (e.g. vinegar or lime scale removers), since they may cause damage to the appliance.
- Keep the appliance away from high temperature and avoid direct sunlight.
- Do not submerge the device in water or any other liquid.

- Wipe the device with a soft cloth and a mild, non-abrasive cleaner when necessary.
- Clean the water tank with soapy water; use a brush to remove any dirt inside and wipe clean with a cloth.
- Ensure the water tank is completely dry before storing the device for long periods.

Precautions

Charging:

- Ensure the device is powered off before connecting the charging cable. Avoid forcefully pulling on the power supply wire or the device.
- Always attach the power cord to the product first, then to the power outlet.
- If the USB cable accidentally falls into water, unplug the adapter first before retrieving the product.
- Do not unplug the adapter with wet hands.
- Do not plug this device into a voltage system that is different from the voltage system specific on the device or charger.
- Prevent overcharging by unplugging the device once the designated charging time is reached.
- Always keep the product away from water sources while charging, and avoid charging it near baths, sinks, or any liquid containers.
- Do not charge on soft surfaces like beds and bedding, especially in enclosed spaces. Should the product overheat during charging, disconnect it immediately.
- Avoid charging in high-temperature areas or near open flames to mitigate fire risks.
- Avoid damaging, breaking, or modifying the USB cable. Do not stretch, bend or twist the USB cable excessively, and do not place heavy objects on it or trap it between objects.
- If the power cord becomes damaged, replace it using only compatible cords or components purchased directly from the original manufacturer.

Usage:

- This product is designed solely for cleaning gum and teeth. Do not use it for any other purpose.
- Do not use the product if you lose the nozzle.
- Only use tips and accessories that are recommended.
- Do not press the nozzle too firmly against your teeth or gums, or it may cause injury.
- People suffering from tooth or gum pain may not be able to use the appliance and should consult their dentist before use.
- Do not insert any foreign object into any opening on the product.
- Do not use the product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and shorten lifespan of the product. Use of these products may void warranty.
- Remove any oral jewellery prior to use of this product.
- Always power off the device prior to cleaning or performing maintenance.
- If the device experiences any malfunction or damage, immediately power it off and avoid unauthorised disassembly.

- This trimmer is waterproof (IPX7 Rating) however do not immerse the product in water to prevent malfunctioning, fire, electric shock or explosion.
- Ensure the product is completely dry before storing to maintain its condition and functionality.
- We recommend avoiding storage in damp, humid areas and high temperatures. Do not throw product into the fire or heat it.

Safety:

- The appliance can be used by children above 8 years of age however supervision is recommended to ensure safe operation. Always store the device out of reach of unsupervised children to prevent any accidental harm.
- Children should not play with the device as this appliance is not a toy.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. Also, the deadly amoeba, naegleria fowleri, may be present in some tap water or non chlorinated well water and may be fatal if directed in these areas.
- Keep the appliance, cable, nozzle out of reach of children under 8 years. The nozzle is small enough to be swallowed and can be lodged in the throat.

Troubleshooting:

Question	Causes	Solutions
Water tank leaks	Seal ring is abnormal	Remove the seal ring and reinstall
Low water pressure	Low battery	Recharge the device
	Nozzle head blocks	Replace with the new nozzle head
	Dirt or debris stuck in nozzle	
Lasts only few minutes even after full charge	Dead battery	Contact service centre for support
	Insufficient charging time	Charge the unit for at least 5 hours
	Charging cord damaged	Try charging the device with new charging cable
	Adapter broken	Try the new adapter
Motor makes sound but water does not come	Water tank is empty	Fill the water tank
	Suction valve dislocation	Keep device on, fill water to the hole where nozzle head goes into

out	Nozzle head blocks	Replace with new tips
Motor doesn't make sound	Low battery	Activate the battery by charging the unit for 16 hours
Product stopped working suddenly	Overheating protection is activated (continuously used 5 tanks of water with highest pressure)	Excessive heat will cause the machine to stop. Switch off the device, cool for 2 hours, and restart

Environmental Protection and Recycling

Safe Appliance Disposal

At the end of its working life, the battery must not be disposed of as urban waste. It must be taken to a special authority differentiated waste collection centre or to a dealer providing this service. Disposing of a household appliance separately avoids possible negative consequences for the environment and health deriving from inappropriate disposal. This enables the constituent materials to be recovered to obtain significant savings in energy and resources.

As mentioned, the built-in rechargeable battery in this product contains substances that could be harmful to the environment. To ensure the responsible disposal and recycling of the device, please follow your country's rules for the separate collection of electrical products and batteries, or refer to the following steps:

1. Activate the device's power switch and let it run until it stops functioning, ensuring the battery is entirely drained.
2. Before disassembling the device and removing its battery, always ensure the device is powered off and disconnected from any power sources.
3. Once the battery has been extracted, never reconnect the device to its charger. Additionally, be cautious to avoid causing a short circuit between the battery's positive and negative terminals.
4. Then follow local authority guidelines to dispose of the removed battery.



Accessories

Need replacement parts? Lost something? Accessories can be purchased via our website www.barberboss.co.uk/shop.

Customer Service

For assistance visit our website: www.barberboss.co.uk or email us at care@barberboss.co.uk. You can also call or WhatsApp us on 0800 612 6688. Or head to our YouTube for some helpful tips/guides www.youtube.com/@BarberBossUK.

