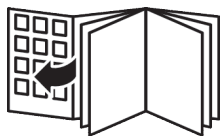
**PHILIPS**

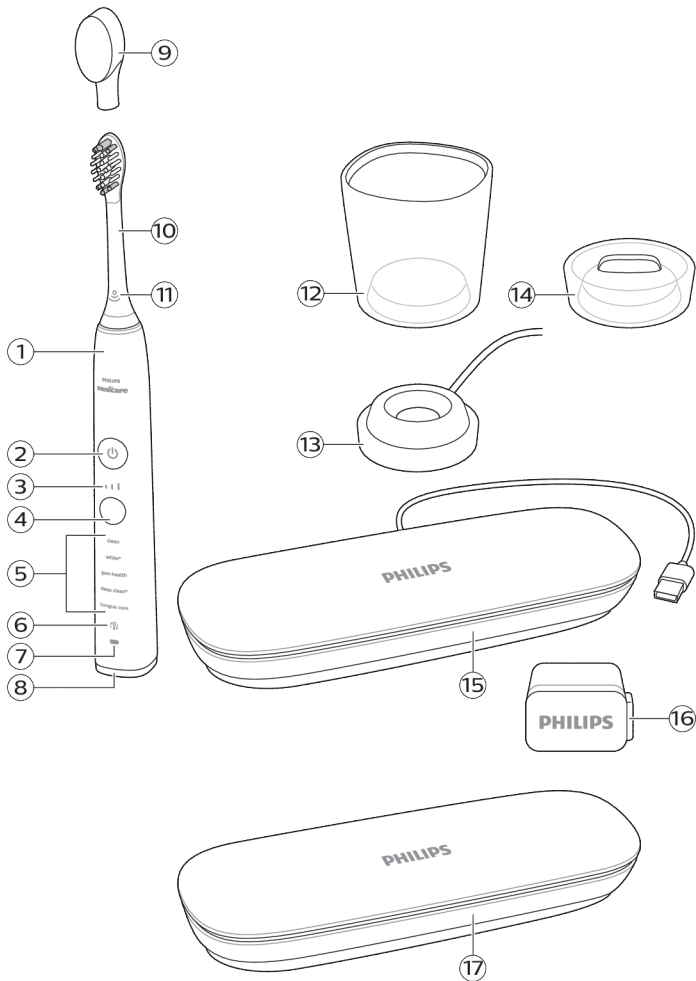
**sonicare**

**9400**

*DiamondClean  
Smart*

English 6  
Azərbaycanca 27  
Čeština 49  
Dansk 70  
Deutsch 91  
Eesti 114  
Español 136  
Français 160  
Hrvatski 183  
Italiano 204  
Latviešu 227  
Lietuviškai 249  
Magyar 271  
Nederlands 293  
Norsk 315  
O'zbek 336  
Polski 359  
Português 383  
Shqip 405  
Slovenščina 427  
Slovensky 448  
Srpski 470  
Suomi 491  
Svenska 513  
Türkmençe 534  
Ελληνικά 555  
Български 580  
Кыргызча 604  
Македонски 627  
Русский 651  
Тоҷикӣ 676  
Українська 700  
Қазақша 724

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## Introduction

Welcome to the Philips Sonicare family!

This toothbrush enables you to obtain superior plaque removal, whiter teeth and healthier gums. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time.

Further support and product registration are available to you at: **[www.philips.com/support](http://www.philips.com/support)**

## Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

## Warnings

- Keep the chargers away from water.
- After cleaning, make sure the chargers are completely dry before you connect USB Charger.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Only use original Philips accessories or consumables. Only use the charger provided with the product. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country. If the cord is damaged, replace the charger with a new one.
- Do not use the chargers outdoors or near heated surfaces.

- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

## **Battery safety instructions**

- This appliance contains batteries that are only to be removed by skilled persons
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables. Only use detachable supply unit CBAyyyy, TCAxxxx, or WAAxxxx.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.

- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on recycling.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
- When you handle batteries, make sure that your hands, the product and the batteries are dry.
- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings). Do not wrap batteries in aluminium foil. Tape battery terminals or put batteries in a plastic bag before you discard them.
- Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).

## **Storage and transportation conditions**

Temperature -10°C to 60°C (14 °F to 140 °F).

## **Electromagnetic fields (EMF)**

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## **Radio Equipment Directive**

Hereby, Philips declares that electric toothbrushes with radio equipment type (Bluetooth, NFC) are in compliance with Directive 2014/53/EU.

- The Bluetooth radio frequency interface in applicable products operates at 2.4 GHz.

- The maximum output power by the Bluetooth appliance is 3 dBm.
- The NFC radio frequency interface in applicable products operates at 13.56 MHz.
- The maximum RF power transmitted by the appliance is 30.16 dBm.

The full text of the EU declaration of conformity is available at the following internet address:

**[www.philips.com/support](http://www.philips.com/support)**

Note: The Features of different models may vary. Some models may not be equipped with Bluetooth or NFC.

## Intended use

Sonicare DiamondClean Smart power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare DiamondClean Smart power toothbrushes are intended for personal use. Use by children should be with adult supervision.

## Your Philips Sonicare (Fig 1.)

- 1 Handle
- 2 Power on/off button
- 3 Intensity indicator
- 4 Mode/intensity button
- 5 Mode indicator
- 6 Brush head replacement reminder
- 7 Battery indicator
- 8 Brushing feedback light

### Accessories:

- 9 Brush head cap
- 10 Smart brush head(s)
- 11 BrushSync symbol
- 12 Charging glass
- 13 Charging base
- 14 Charging stand
- 15 Charging travel case with integrated USB cord
- 16 USB-A wall adapter
- 17 Travel case

Note: Accessories included may vary based on the model purchased.

## Sonicare app - Getting started

The Sonicare app pairs with your toothbrush to provide you a connected experience. By connecting your toothbrush to your Sonicare app account, you will be able to:

- Track your brushing progress.
- Receive personalized tips and actionable recommendations to improve your oral health.
- Access the full range of benefits and receive ongoing upgrades.

The app is compatible with a wide range of iPhone and Android™ smartphones.

To start using the Sonicare app:

- 1 Download the Sonicare app to your phone.
- 2 Ensure your phone's Bluetooth is turned on.
- 3 Pick up your toothbrush to ensure it is active (lights on).
- 4 Open the app and follow the steps.
- 5 Pair your toothbrush with the app.
- 6 Create your account via the app. Complete firmware update, if prompted, to access the latest improvements and features.
- 7 Brush regularly with your toothbrush connected to the app. When you regularly sync the toothbrush with the app, you will receive updates to help improve your oral healthcare.
- 8 Sync regularly.
- 9 Keep your phone close to your toothbrush for optimal data transfer.

**To sync manually:** Pair/connect your toothbrush with the app every two weeks to benefit from the app features/capabilities.

**To sync automatically:** Allow location permission while setting up the app. By allowing location permission, your phone knows when it is in the connection range of your toothbrush and can refresh your brushing data in the app to provide the latest insights and recommendations.

Note: Make sure your phone's Bluetooth is switched on when using the app so that your toothbrush can transfer and update your brushing data in the app. If you have questions about why your brushing data is collected, be sure to review the Sonicare Privacy Statement, available throughout the app setup process.

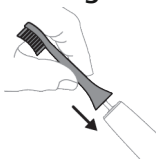
## App features

Your Sonicare power toothbrush is equipped with smart sensors to provide you with feedback about your brushing habits, including:

- Daily brushing sessions
- Time spent during each brushing session
- Applied brushing pressure
- Personalized tips and actionable recommendations
- Brush head replacement reminder (based on actual use)

The Sonicare app collects the data provided by the smart sensors and tracks your brushing habits over time.

## Getting started

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.

Note: It is normal to see a slight gap between the brush head and the handle.

## Using your Philips Sonicare toothbrush

If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first-time users to apply too much pressure. Only apply gentle

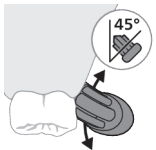
pressure and let the toothbrush do the brushing for you. Please follow the brushing instructions below for the best experience.

To help you transition to your new Sonicare power toothbrush, use a lower intensity setting for the first 1-2 weeks then transition to a higher setting.

## Brushing instructions


- 1 Wet the bristles and apply a small amount of toothpaste.


- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees). Pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

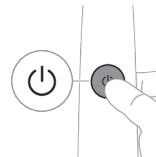Note: Keep the center of the brush in contact with the teeth at all times.

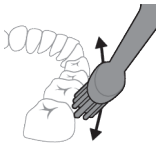
- 3 Press the power on/off button to switch on the Philips Sonicare.

- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

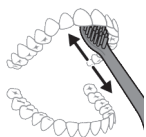
Note: If you apply too much pressure, you will be alerted by a change in vibration of the handle and a flashing purple brushing feedback light.

- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.
- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush switched on or off, as you prefer.

Note: Your Philips Sonicare toothbrush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised.

## BrushPacer

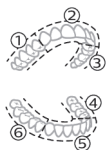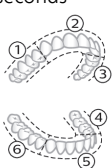
To make sure you brush evenly throughout the mouth, divide your mouth into 6 sections using the BrushPacer feature. The BrushPacer divides the brushing time into 6 equal segments and indicates when you should move to the next area. The segments are indicated with a brief pause in vibration. The toothbrush automatically stops at the end of the brushing session.

Note: Additional segments may be added when you use the White+ or Gum Health Modes.

## Modes

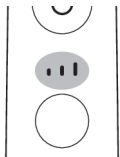
Your power toothbrush is equipped with up to 5 different modes: Clean, White+, Gum Health, Deep Clean+ and Tongue Care (specific models only). Based on the smart brush head you use, the optimal mode and intensity will be selected automatically by the BrushSync technology. To change modes, press the mode/intensity button while the handle is turned off.

	<b>Clean and Deep Clean+</b>	<b>White+</b>	<b>Gum health</b>	<b>Tongue Care (specific models only)</b>
<b>Benefit</b>	Plaque removal	Plaque and surface stains removal	Plaque removal and gum massage	Tongue cleaning

	<b>Clean and Deep Clean+</b>	<b>White+</b>	<b>Gum health</b>	<b>Tongue Care (specific models only)</b>
<b>Recommended intensity level</b>	High	High	High	Low
<b>Total brushing time</b>	Clean: 2 minutes Deep Clean+: 3 minutes*	2 minutes and 40 seconds	3 minutes and 20 seconds	20 seconds
<b>How to brush</b>	<p>Clean: Brush each segment for 20 seconds</p>  <p>Deep Clean+: Brush each segment for 30 seconds*</p> 	<p>Brush each segment for 20 seconds,</p>   <p>then brush segments 7 and 8 for 20 seconds each.</p>  	<p>Brush each segment for 20 seconds,</p>   <p>then brush segments 7, 8, 9 and 10 for 20 seconds each.</p>  	<p>1. Wet the micro bristles and brush your tongue with a back and forth motion. Clean as far back as comfortable. 2. Rinse your mouth and repeat the procedure two more times.</p>

\*Note: When brushing with the Philips Sonicare app in 'Deep Clean+' mode, the brushing time will be 2 minutes when connected to the app and 3 minutes when not connected to the app.

## Intensities

Your power toothbrush comes with 3 different intensity settings:

- low intensity (one bar)
- medium intensity (two bars)
- high intensity (three bars)

When attached to the handle, the smart brush head will automatically select the recommended intensity. If you wish to change to a different intensity, you can do so by pressing the mode/intensity button while brushing. The intensity cannot be changed as long as the handle is powered off or paused.

## BrushSync technology

BrushSync technology enables your brush head to communicate with your handle using a microchip. The  symbol at the bottom of the brush head indicates that the brush head is equipped with this technology.

BrushSync technology enables:

- Brush head replacement reminder
- BrushSync mode pairing

Philips Sonicare offers a wide range of smart brush heads, equipped with BrushSync technology. To explore our full range of toothbrush brush heads, and find the brush head that is best for you, go to

**[www.philips.com/toothbrush-heads](http://www.philips.com/toothbrush-heads)** for more information.

## Features

- Pressure Sensor Feedback
- BrushPacer
- SmartTimer
- Brush Head Mode Pairing
- Easy Start
- Brush Head Replacement Reminder

## Pressure Sensor Feedback

Your Sonicare toothbrush measures the pressure you apply while brushing, to protect your gums and teeth from damage.

If you apply excess pressure, the handle changes its vibration and the brushing feedback light flashes purple at the bottom of the handle until you reduce the pressure.

**Note:** Your product comes with the Pressure Sensor activated. To deactivate this feature (see 'Activating or deactivating features').

## BrushPacer

The BrushPacer is an interval timer that helps you to brush evenly through the mouth by dividing your mouth into 6 sections. (see 'Brushing instructions')

## SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

**Note:** If you press the power on/off button after you started the brushing cycle, the toothbrush pauses. After a pause of 30 seconds, the SmarTimer resets.

## Brush Head Mode Pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates with the handle and automatically pairs it with the recommended mode and intensity.

Most brush heads with BrushSync technology will auto pair to Clean Mode. Brush heads with specific benefits will pair with the corresponding mode (when available), examples include:

- G3 Premium Gum Care : Gum Health mode
- W3 Premium White: White+ mode
- TongueCare+: Tongue Care mode

**Note:** If you choose to change mode/intensity from the recommended setting, the brush head will remember the your chosen setting for future brushing sessions.

## EasyStart

This Philips Sonicare model comes with the EasyStart feature deactivated. To activate EasyStart (see 'Activating or deactivating features'). The EasyStart feature gently increases the power over the first 14 brushings to help you get used to brushing with the Philips Sonicare. The EasyStart feature works with all modes.

### Disclaimer

When the Philips Sonicare toothbrush is used in clinical studies, it must be used in clean mode at high intensity with the handle fully charged and the EasyStart feature deactivated.

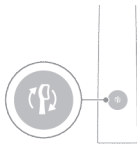
## Brush Head Replacement Reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

- 1 When you attach a new smart brush head for the first time, the handle recognizes that you have a Philips brush head with BrushSync technology and begins tracking brush head wear.
- 2 Over time, based on the pressure you apply and the total time the brush head has been used, the handle tracks brush head wear in order to determine the optimal time to replace your brush head. This feature guarantees the best cleaning and care of your teeth.
- 3 When the brush head replacement reminder indicator lights up amber, you should replace your brush head.

All smart brush heads come with the brush head replacement reminder feature activated. To deactivate the brush head replacement reminder see 'Activating or deactivating features'.


## Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush Head Replacement Reminder
- Pressure Sensor Feedback

Note: You can only activate or deactivate the Brush Head Replacement Reminder when a smart brush head is on the handle.

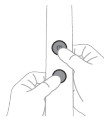Step 1: Place the handle on the charging stand or in the charging glass.

Step 2: Press and hold the mode/intensity button.

EasyStart

Pressure Sensor  
Feedback

Brush Head  
Replacement Reminder

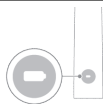
Press the power  
on/off button 1  
time.

Press the power  
on/off button 2  
times.

Press the power on/off  
button 3 times.


Step 3: Let go of the mode/intensity button.

If you see the battery indicator flash white 2 times and hear 2 beeps then the feature has been activated.

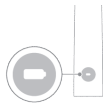OR

If you see the battery indicator flash amber 1 time and hear 1 beep, then the feature has been deactivated.

## Charging and battery status

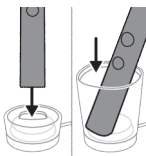
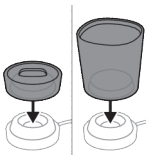
This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day) on a fully charged battery. The battery indicator shows the battery status when you complete the 2-minute brushing session, when you pause the toothbrush, or while it is charging.

Note: It can take up to 24 hours to fully charge the battery before first use, but you can use the Philips Sonicare toothbrush before it is fully charged.

## Charging

### Charging on charging base

- 1 Plug the cord of the charging base into an electrical outlet.
- 2 Place the charging stand or the charging glass on the charging base.
- 3 Place the handle on the charging stand or in the charging glass.



- The battery indicator and light ring illuminate white and the handle beeps twice. This indicates the toothbrush has started charging.
- 4 The charger will charge your toothbrush until it is fully charged.
    - The handle is charging when the battery indicator is flashing white. When the handle is fully charged, the battery indicator lights up white for 30 seconds and then switches off.

### Charging with the travel case

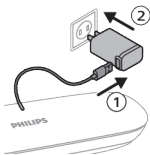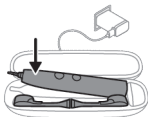Note: Accessories included may vary based on the model purchased.

- 1 Take the USB cable from the bottom of the charging travel case.
- 2 Connect the charging travel case via the USB cable to 5V DC input power either from:
  - The USB adapter provided
  - a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
  - a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
  - a UL-listed or certified 1310 Class 2 power supply.

- 3 Place the handle in the travel case. The travel case charges your toothbrush until your toothbrush is fully charged.

## Battery status

### Battery status (handle on plugged in charger/charging travel case)

#### Handle on plugged in charger/charging travel case:

When the handle is placed on the charger or in the charging travel case, the battery indicator shows the battery level.

Battery status	Battery indicator
Charging	Flashes white
Full	Solid white for 30 seconds and then switches off

### Battery status (when handle is not placed on charger/charging travel case)

#### Handle is not placed on charger/charging travel case:

When the toothbrush is awake, the battery indicator at the bottom of the handle shows the battery level.

Battery status	Battery indicator	Sound
Full	Solid white	-
Partially full	Flashes white	-
Low	Flashes amber	The handle emits 3 beeps after brushing
Empty	Flashes amber	The handle emits 2 sets of 5 beeps after brushing

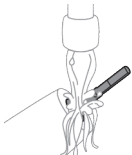
Note: To save energy, the battery indicator will switch off when not in use.

Note: When the battery is completely empty, the toothbrush switches off. Place the toothbrush on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your toothbrush on the charger when not in use.

## Cleaning

### Toothbrush handle


- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

**Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.**

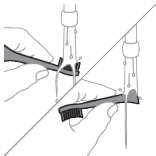- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

### Brush head

- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

### Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

### Travel case

Rinse with warm water and use a damp cloth to clean inside of travel case.

### Charging travel case

- 1 Unplug the charging travel case before cleaning.
- 2 Use damp cloth to wipe the surface of the travel case.

## Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

## Replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

## Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX99XX).

## Warranty and support

If you need information or support, please visit **[www.philips.com/support](http://www.philips.com/support)** or read the international warranty leaflet.

## WARRANTY EXCLUSIONS

What is not covered under warranty:

- In mouth attachments including brush heads and nozzles.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

## Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

The built-in rechargeable battery **must only be removed by a qualified professional** when the appliance is

discarded. Before removing the battery, make sure that the battery is completely empty.

## Removing the rechargeable battery

The built-in rechargeable battery **must only be removed by a qualified professional** when the appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.**

**When you handle batteries, make sure that your hands, the product and the batteries are dry.**

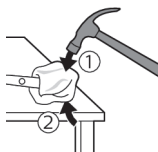
**To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.



- 3** Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

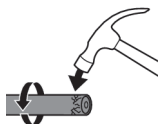
Note: You may have to hit on the end several times to break the internal snap connections.

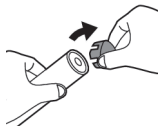
- 4** Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

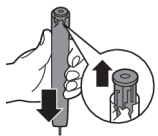
- 5** Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

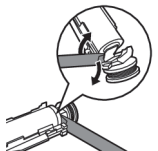- 6** Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

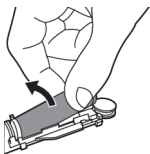- 7** Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.


- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

**Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.**

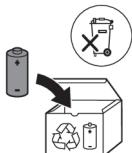- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

**Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.**

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©2022 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.