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
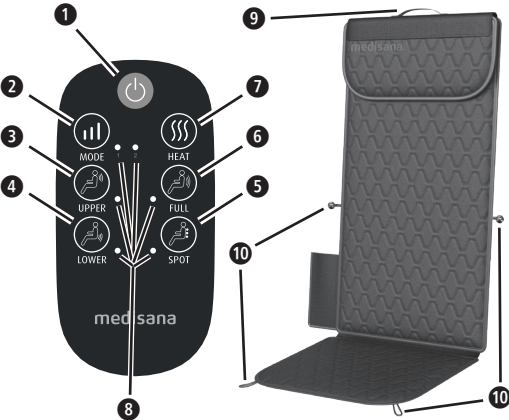
Device and controls

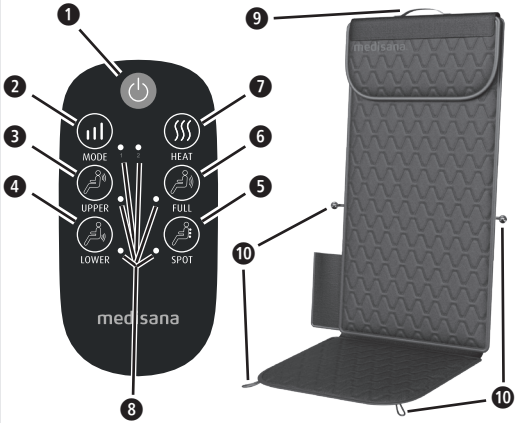Explanation of symbols

The device must only be used in enclosed spaces!

Do not use the device in damp rooms (e.g. when bathing or showering).

IMPORTANT
Non-observance of these instructions can result in serious injury or damage to the device.

WARNING
These warning notes must be observed to prevent any injury to the user.

CAUTION
These notes must be observed to prevent any damage to the device.

NOTE
These notes give you useful additional information on the installation or operation.

Protection category II

LOT number

Manufacturer

Recycling symbols / codes:
These are used to provide information about the material and its proper use and recycling.

Do not wash!

GB Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- The device must only be used for its intended purpose as described in the instruction manual. Using the appliance for any other purpose invalidates the warranty.
- The device is not to be used for commercial use or medical purposes.

Safety Instructions – About the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only insert the power supply unit into a socket when the device is switched off.
- This device may only be operated with the mains adapter (XZ1200-2500G) supplied.
- Keep the power supply unit, cable and appliance well away from heat, hot surfaces, moisture and liquids.
- Do not touch the power supply unit if you are standing in water or if your hands are wet or damp.
- Do not take hold of a device if it falls into water. Pull out the power plug and the power supply unit from the socket immediately.
- Only plug the device into a socket where the plug is freely accessible.
- Always switch off the appliance after use with the buttons on the control unit and unplug the power supply unit.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the device by the cable or the power supply unit.
- You must not continue to use the unit if the cable or power supply unit are damaged. For safety reasons, these parts may only be replaced by authorized service personnel. In order to prevent hazards, always send the unit in to the service centre.
- Ensure that no one can trip over the cable. The cable must not be kinked, trapped or twisted.

Safety Instructions – Specific groups of people

- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Do not use this device to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the device, or only after consulting your doctor, if:
 - » you are pregnant,
 - » you have a cardiac pacemaker, artificial joints or electronic implants,
 - » you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Take care when using the heating function. The device has hot surfaces. Do not use the device on infants, vulnerable individuals or people who are not sensitive to heat.
- Before using the device, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

Safety Instructions – Before using the unit

- Check the power supply unit, cable and Massage Seat Cover carefully for damage before each use. Do not attempt to use a defective unit.
- Do not use the device if it or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- Do not keep or place any electrical device on a surface where it is in danger of falling into a bathtub or sink.
- Only use the device in the proper place, as described in this instruction manual.
- Do not use the device while you are asleep or when you are lying in bed.
- Do not use the device before going to sleep. The massage has a stimulating effect.
- Never use the device while driving a vehicle or operating a machine.
- The maximum operating time for one session is 15 minutes. Longer use shortens the life of the device and can lead to tension of the muscles instead of relaxation.
- Very long periods of use can lead to overheating. Always allow the device to cool down before using it again.
- Never store or use the device next to an electric stove or other sources of heat.
- Never leave the device unattended if it is connected to the power supply.
- Do not allow the device to come into contact with pointed or sharp objects.
- Do not cover the device when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.
- Never sit or stand on the moving parts of the device with all your weight and do not place objects on the device.
- Do not use an external time switch or a separate remote control system to operate the device.

- Do not use the device whilst using aerosol sprays or simultaneously using oxygen.

Safety Instructions – Maintenance and cleaning

- The device is maintenance free. Do not try to open it. There are no parts inside that are to be cleaned or serviced by the user.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all warranty rights. Consult your specialist dealer and have repairs carried out by the service centre to avoid dangers.
- Do not immerse the device in water or other fluids. If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.
- Cleaning and user maintenance shall not be carried out by children without supervision.

Intended use

The device is used to massage certain body parts for relaxation and well being. It is not intended for medical therapeutic use. It is not suitable for clinics, massage practices, wellness studios etc. The massage device is not a medical device and is not intended for the treatment of pain. In case of doubt or uncertainty, consult your doctor before use.

Scope of delivery

Please check first whether the unit is complete and does not show any damage. If in doubt, do not operate the unit and contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Massage Seat Cover **MC 700**
- 1 Mains adapter XZ1200-2500G (cable length: approx. 180 cm)
- 1 Quick Start Guide (short instructions)

If you notice any transport damage when unpacking, please contact your dealer immediately.

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.

WARNING
Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Device and controls

- ➊ On/Off button
- ➋ Mode button (slow = 1 / fast = 2), indicated via LEDs “1” or “2”, resp. Massage type buttons (indication via the LED positioned next to the respective button):
 - ➀ upper back
 - ➁ lower back
 - ➂ spot massage
 - ➃ full back
- ➌ Heat button (If active, the button lights up in red.)
- ➍ LEDs
- ➎ Carrying handle
- ➏ The retaining loops and pearls secure the seat cover when folded.

Using the device

- Place the Massage Seat Cover on a high chair or armchair. The handle ➉ allows for easily carrying the device. Ensure that the device is supported over its entire area by a sufficiently high backrest.
- Alternatively, you can lay the **MC 700** flat on a mattress or couch (so that you can recline on it).
- Insert the plug into the mains socket. Ensure that the power supply is easily accessible.
- Take a seat on the chair, sit comfortably and upright (or lay down on the mattress / couch). Ensure your position is centered and such that the massage heads can move to the left and right of your spinal column.
- Familiarise yourself with the functions.
- Use the buttons on the integrated controller panel to operate the device. Switch it on by pressing the On/Off button ➊. The massage starts in operating mode No. 1 (“slow”) and “full back”.
- A short press on the Mode button ➋ or any massage type button ➀...➃ switches to the respective operating mode. The active mode is indicated by the LEDs ➍.
- Press the Heat button ➌ to switch on the heat function. The button then lights up in red. The heat function can only be used during a massage. Pressing the button ➌ again deactivates the heat function.
- Pressing the On/Off button ➊ again switches the device off.
- If the On/Off button ➊ is not pressed, the device switches off automatically after 15 minutes.
- When the device switches off (either automatically or manually by pressing the On/Off button ➊) the massage heads are moved downwards into their parking position.
- Always allow the appliance to cool down completely before using it again.

Cleaning and maintenance

- Before cleaning the unit, make sure it is switched off and disconnected from the mains. Clean the unit once it has cooled down using a soft, damp cloth.
- Under no circumstances should you use aggressive detergents, strong brushes, solvents or alcohol, as they may damage the surface of the unit.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- You should ideally store the device in its original packaging and keep it in a clean, cool and dry place.

Technical Data

Name and model:	medisana Massage Seat Cover MC 700
Power supply:	via mains adapter, type XZ1200-2500G, input: 100–240V~ 50/60Hz, output: 12V≐ 2.5A 30W, when being switched off: <0.3W 15 minutes, with automatic switch-off approx. 78 × 47.4 × 6.5 cm approx. 2.5 Use only in dry indoor areas Cool, clean and dry 88968 40 15588 88968 4
Power consumption:	
Max. operating time:	
Dimensions (L×W×H):	
Weight:	
Operating conditions:	
Storage conditions:	
Item number:	
EAN number:	

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:


1. The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH, Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under <https://docs.medisana.com/88968>.

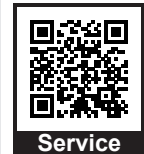Information about service can be found here:
www.medisana.com/servicepartners

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