emmi[®]DENT

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emmi[®]-dent SLIM

Emmi Ultrasonic GmbH
Gerauer Strasse 34
D-64546 Mörfelden-Walldorf
Tel: +49 (0)6105 - 40 67 94
Email: info@emmi-dent.de

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Warranty services

Warranty period:

Your ultrasonic toothbrush is guaranteed for 2 years from the date of purchase.

Scope of guarantee:

This warranty covers material and manufacturing defects on the device. It includes the motor and the battery.

Not covered:

Normal wear and tear, damage caused by improper use, accidents, modifications or unauthorised repairs are excluded from the warranty.

Warranty services:

In the event of a warranty claim, we offer to repair or replace the appliance.

Prerequisites for utilisation:

To make a warranty claim, please keep the proof of purchase and contact our customer service.

Rights of the consumer:

This guarantee does not limit your statutory rights, but supplements them.

© Copyright emmi EMAG AG. All rights reserved. Manufactured in Germany/EU. The contents of these operating instructions are subject to copyright. The operating instructions may not be copied in whole or in part. Modifications, abridgements, extensions and supplements as well as any publication or translation by third parties require the prior written consent of emmi EMAG AG.

emmi® ULTRASONIC

Emmi Ultrasonic GmbH • Gerauer Strasse 34

64546 Moerfelden-Walldorf • Germany

Tel. +49 (0)6105 40 67 94 • Fax +49 (0)6105 40 67 81

www.emmi-dent.com • info@emmi-dent.de

Tips for optimal dental care

Use dental floss (emmi® DENT dental floss) and mouthwash (emmi® DENT mouthwash) in addition to the ultrasonic toothbrush. Visit your dentist regularly for check-ups and professional dental cleanings.

The brush head has the following features:

- Previous mode is saved (mode 1, 2 or 3)
- When to change the ultrasonic attachment?
Your intelligent attachment for the ultrasonic toothbrush memorises the entire period of use and reminds you in good time when it needs to be replaced. This ensures optimum dental care at all times. Look out for the LED indicator with the attachment, which tells you to change the attachment.

- Stabilises the ultrasonic power
Active communication between the intelligent attachment and the handpiece of your ultrasonic toothbrush ensures stabilisation of the fluctuating ultrasonic pulses. This advanced technology ensures even and effective cleaning of your teeth. Make sure that the intelligent attachment is correctly positioned on the handpiece to ensure the optimum performance and efficiency of your toothbrush.

Thank you for choosing our ultrasonic toothbrush.
We wish you a thorough and pleasant dental cleaning!

3. brushing teeth

Switch on the ultrasonic toothbrush by pressing the on/off switch. Select the desired mode with the setting button. The intelligent brush head memorises the selected mode. Hold the toothbrush on your teeth and move it from tooth to tooth every 6 seconds. You will hear a haptic signal every 6 seconds to help you change the brushing position. Brush your teeth for at least the three minutes recommended by dentists. The appliance switches off automatically after 6 minutes. This is the travel protection so that the battery is not discharged unintentionally. If no attachment is attached, the appliance switches off within one second and you cannot brush your teeth.

4. after use

Switch off the ultrasonic toothbrush. Rinse the brush head thoroughly under running water. Place the ultrasonic toothbrush on the charging station when the battery indicator (red light) signals that it needs to be charged.

5. care and maintenance

Cleaning:

Clean the handle of the toothbrush regularly with a damp cloth. Brush heads should be replaced every three months or when wear is visible and noticeable (the lack of smoothness on the teeth).

Storage:

Store the toothbrush in a dry place when not in use.

3. adjust settings

The toothbrush has three different modes that can be selected using the setting button (top button).

Mode 1:

Light vibration for sensitive teeth and gums.

Mode 2:

Medium vibration for healthy teeth and gums.

Mode 3:

Strong vibration for increased ultrasonic power and fast cleaning.

Press the setting button to switch between modes. Each mode is signalled by a different LED display:

Mode 1:

Yellow light

Mode 2:

Purple light

Mode 3:

Turquoise light

Instructions for use

1. attach the brush head

Attach the brush head firmly to the handle of the toothbrush (handpiece).

2. apply toothpaste

Apply a pea-sized amount of ultrasonic toothpaste (emmi DENT ultrasonic toothpastes) to the bristles.

emmi®-dent - The original for innovative tooth cleaning with 100% original ultrasonic technology

Welcome and congratulations on your new generation ultrasonic toothbrush with intelligent brush head! These instructions will help you get the most out of your ultrasonic toothbrush.

Scope of delivery

- 1x ultrasonic toothbrush (handpiece)
- 1x charging station + mains adapter (USB-C)
- 2x intelligent brush attachments
- 1x ultrasonic toothpaste

First steps

1. unpacking and charging


Carefully remove the ultrasonic toothbrush and accessories from the packaging. Place the ultrasonic toothbrush on the charging station. Connect the charging station to a power socket and allow the toothbrush to charge fully (approx. 24 hours before first use).

Red light: The ultrasonic toothbrush is charging.

Green light: The ultrasonic toothbrush is fully charged.

2. switch on/off

Press the on/off switch (bottom button) once for approx. 0.5 seconds to switch on the ultrasonic toothbrush. Press the same button again to switch the toothbrush off.

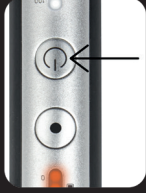
emmi! DENT

INSTRUCTION MANUAL

emmi![®]-dent SLIM

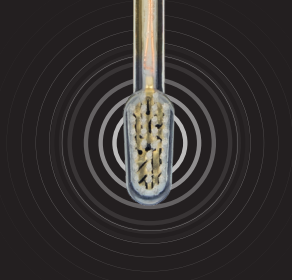with **100%**
ULTRASONIC

Emmi Ultrasonic GmbH
Gerauer Strasse 34
64546 Moerfelden-Walldorf
Germany
Tel: +49 (0)6105 - 40 67 94
Email: info@emmi-dent.de