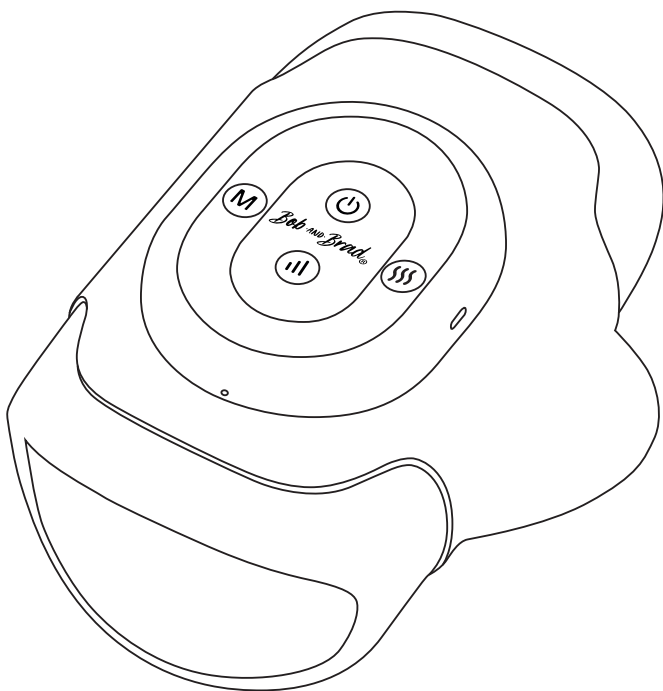
User Manual

HandSpa Pro

Hand Massager


Model: H50

Meet Bob and Brad

Bob and Brad are dedicated to creating a healthier world. We bring to you some amazing products that we have personally vetted and truly believe in! We want only the best for you and the best is right here! Enjoy!

NEED HELP?

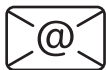SUPPORT@BOBANDBRAD.COM

www.bobandbrad.com

Tel: +1(612) 567-3035

M-F 9:00 AM- 5:30 PM

Your guide to being healthy, fit, and pain-free

B&B Innovation LLC
3310 S Fairview St, Santa Ana, CA 92704

Maintenance

1. Turn off the device before cleaning it. Wipe the device with a soft, damp cloth and non-abrasive cleaning agents.
2. Do not immerse the device in liquid.
3. Do not expose the device to direct sunlight or high temperatures. If the device has not been used for over 3 months, recharge it before use.
4. Do not disassemble the device yourself. Otherwise, it may cause damage and void the manufacturer's warranty.
5. Battery life may vary depending on usage and operating conditions. To prolong its lifespan, if the device will not be used for an extended period, fully charge the battery and keep the device powered off. Even when not in use, recharge the device every 2 months.
6. If the product is not used for an extended period, store it in a PE bag and keep it in a dust-free environment.

Warnings and Safety Cautions

Please read this user manual carefully and use the product according to the instructions.

DANGER - To reduce the risk of electric shock:

1. Always unplug the charger from the electrical outlet after charging the device.
2. Never use pins or other metallic fasteners on this product.

WARNING - User Precautions:

1. If you have any of the following conditions, consult your doctor before using this product:

- Use of a pacemaker or similar implanted medical device;
- Heart disease, stroke, or high blood pressure;
- Hand disabilities or a history of hand fractures;
- Weak or damaged bones, especially in the hands;
- Recent surgery or surgeries;
- Pregnancy, recent childbirth, or menstruation;
- Malignant tumors;
- Osteoporosis;
- Blood circulation problems or blood clots;
- Skin infections, skin diseases, or bruises;
- Ongoing medical treatment or feeling unwell;
- Limb abnormalities;
- Previous hand injuries or ailments due to illness or accidents;
- Diabetes mellitus, severe peripheral circulatory disorders, sensory impairment, or anemia.
- If you feel unwell or experience any unusual discomfort, stop using the device immediately.

2. Children, people with disabilities, and the elderly should only use this product under proper supervision.
3. This device includes a heating function. Users with reduced sensitivity to heat should exercise extra caution.
4. This product is intended for personal use only and is not a substitute for medical care or medical devices.
5. Do not wear rings, watches, bracelets, jewelry, or other hand accessories while using the massager.
6. Do not use the device for more than 15 minutes at a time.
7. Do not disassemble the device.
8. Do not expose the lithium battery in the device to heavy impact, high pressure, or fire.
9. Do not use the device if its battery has been exposed to water.
10. Do not use the device in high-temperature or humid environments.
11. Only use a DC 5V 1A or DC 5V 2A charger (not included) to charge the device.
12. Keep the device away from sharp objects.
13. Do not use the device while driving or walking.
14. Do not use the device while it is charging.
15. Do not use the device with wet hands.
16. Do not use the device if it is not functioning properly, has been damaged, or has been submerged in water.
17. Do not use the device if, during battery use or charging, it emits an odor, generates excessive heat, or shows signs of deformation.
18. Do not use the device under a blanket or pillow, as this may cause overheating, fire, electric shock, or injury.

FCC Statement

This device complies with Part-15 of the FCC Rules.

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference.
- (2) This device must accept any interference received, including interference that may cause undesired operation.

FCC Statement NOTE

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause any harmful interference to radio or television receptions, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help with the supplied fastener.
- DO NOT hang the unit by the power cord.

Disclaimer

Please contact your physician before using the Bob and Brad Massager if you have any physical health concerns or have any of the following conditions:

- People with implanted medical electronic devices that are easily disturbed, such as cardiac pacemakers.
- People with malignant tumors.
- People with heart disease.
- People with impaired temperature sensation.
- Women with severe peripheral circulatory nerve disorders caused by diabetes or those in the immediate postpartum period.
- People with sensory disorders.
- People experiencing acute pain due to conditions such as osteoporosis, sprains, or muscle strains, who require rest.
- People with a body temperature exceeding 38°C, such as those with severe symptoms of acute inflammation (e.g., fatigue, chills, blood pressure fluctuations) or those with weakened physical strength undergoing treatment at other medical institutions.
- People with tenosynovitis.


THIS PRODUCT IS NOT A MEDICAL DEVICE.

IC Warning

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:


- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

Main Structure

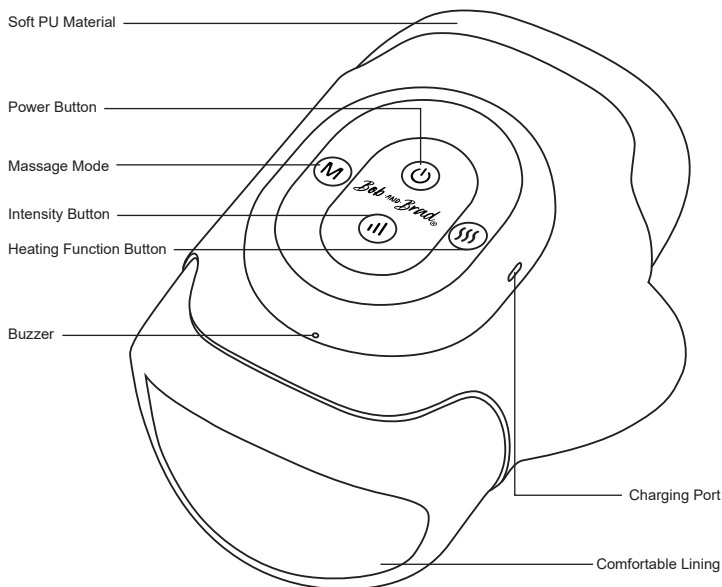
Caution Before Using


This device is suitable for both left and right hands.

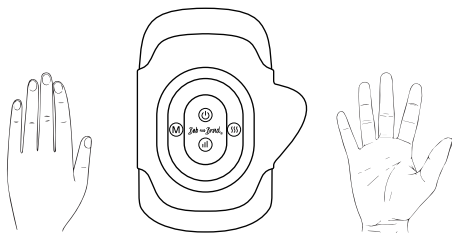
You can use disposable gloves and apply hand cream. Disposable gloves are provided in the package.

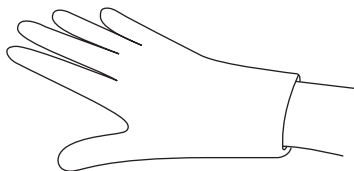
1. Please ensure that the product has sufficient power. If the battery is low, the device will automatically shut down. If it cannot be turned on, charge it before using it again.
2. Remove all accessories on your hand.

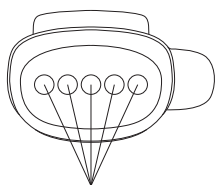Insert your fingers into the corresponding holes.

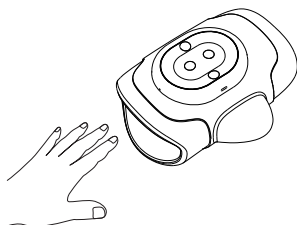Place your hand into the hand massager, ensuring each finger is deeply inserted into the designated finger slots.

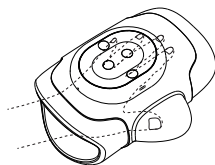Always perform the actions of inserting and removing your hand when the device is powered off.

Button Function Description

Power Button

Press and hold for 1.5 seconds to turn the device on/off.

Low Battery Indicator: The red light flashes.

Charging Indicator: The red or orange light breathes.

Full Charge Indicator: The button light turns green.

Mode Button

Short press to switch the massage mode:

Mode 1 (Vitality Mode + Low Level Intensity): Cyan light

Mode 2 (Relaxation Mode + Medium Level Intensity): Green light

Mode 3 (Intensive Mode + High Level Intensity): Red light

Intensity Adjustment Button

Short press to cycle through intensity levels:

Low Level: Cyan light

Medium Level: Green light

High Level: Red light

Heating Indicator Light


Short press to turn on the heating function. The red light stays on, indicating a constant temperature of 102.2-113°F.

Short press again to turn off the heating function. The red light turns off.

Indicator Light Instruction

 Device Power	<9% Power	Red Flashing Light
	9%~50% Power	Orange Steady Light
	51%~100% Power	Green Steady Light
 Charging	<9% Power	Red Breathing Light
	9%~50% Power	Orange Breathing Light
	51%~99% Power	Green Breathing Light
	100% Power	Green Steady Light
M Massage Mode	Vitality Mode	Cyan Steady Light
	Relaxation Mode	Green Steady Light
	Intensive Mode	Red Steady Light
 Intensity Adjustment	Low Level	Cyan Steady Light
	Medium Level	Green Steady Light
	High Level	Red Steady Light
 Heating Function	Turn On Heating	Red Steady Light
	Turn Off Heating	No Light

Charging Instruction

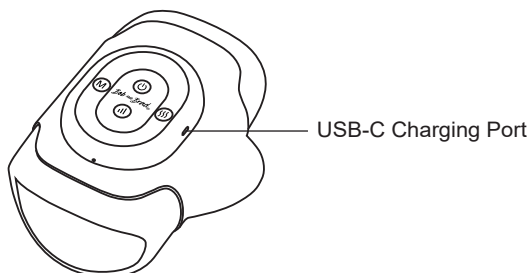1. Use the included charging cable to connect the product to a USB port for charging.
2. The product requires approximately 4 hours to fully charge.
3. When the battery level is below 9%, the power indicator will flash red as a low battery warning. The device will automatically shut down after 1 minute and can only be restarted after charging.
4. The default working time of the product is 15 minutes. To use it again, the device must be restarted.
5. The product cannot be turned on during the charging process.
6. During charging, the indicator light on the power button will display a breathing light effect until it turns solid green when fully charged. Adapter not included.
7. Supports multiple charging methods, including computer USB ports, power banks, and 5V/1A power adapters.

**Use with 5V  1A
charger output**

5V  1A

**Do not use the product
when it is being charged**

Technical Details

- Name: HandSpa Pro Hand Massager
- Model No.: H50
- Input: 5V  1A
- Rated Power: 5W
- Dimensions: 9.3" x7"x 4.4" / 235x178x113mm
- Weight: 728g
- Battery Rated Voltage: 3.7V 
- Battery Capacity: 2500mAh
- Charging Time: 4 Hours

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massager doesn't turn on.	Power is off.	Long press the Power Button to turn on the massager.
	Battery is too low.	Charge the massager.
	Massager shuts-off after 15 minutes.	Long press the Power Button to turn on the massager for another 15 minutes.
Massager shuts off automatically.	The massager's overheat protection will switch off the massager automatically when too hot.	Wait for the massager to cool down before powering it back on.
	Massager shuts off after 15 Minutes.	Long press the Power Button to turn on the massager for another 15 minutes.
The device can't be charged.	The voltage and current do not match the product requirements.	Use with 5V/1A charger output.
	The charging port is not in contact with the power supply.	Make sure the charging port and power supply are properly inserted.
If you encounter any other problems, contact: support@bobandbrad.com		

Guarantee

This product comes with a two-year warranty from the date of purchase. Warranty will be voided due to any of the following:

- Damaged due to improper usage or storage.
- Damaged by attempting to repair or disassemble the device.
- Period of warranty exceeded.

WEEE Instructions

Correct Disposal of this product	
	This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal. Recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.


Battery Instructions

- Do not submerge the product into water or get it wet.
- Do not charge the product near an open flame or in any extremely hot conditions.
- Do not use or store the product near heat sources like furnaces or heaters.
- If the battery leaks or emits an odor, remove it immediately from close proximity to an open flame.
- Do not put the product into a fire or heat it up in any way.
- Do not disassemble the battery in any way.
- Do not use the product if, while the battery is in use or charging, it emits an odor; gives off excessive heat; shows signs of deformation, discoloration, or any other abnormalities. If the battery is in use or being charged, disconnect the battery from the power source.
- Charge for about three hours before its first use to a full charge. The device can be connected or disconnected.
- Connect the output end of the connection cable to the bottom of the massage instrument's DC charging port. Then, connect the input end of the connection cable to a mobile power adapter or cell phone charger (5V/1A; if less than 1A, the charging time will be longer).
- The power indicator light will flash regularly, indicating that the battery is being charged.

RoHs Declaration

Our Product is compliant with RoHS 2.0 (2011/65/EU Directive) regarding restriction of hazardous substances, the raw materials and components of our product meet the following toxic and hazardous substance limits.

Cadmium (Cd): < 100 ppm	Lead (Pb): < 1000 ppm
Mercury (Hg): < 1000 ppm	Hexavalent Chromium: (Cr VI) < 1000 ppm
Polybrominated Biphenyls (PBB): < 1000 ppm	Polybrominated Diphenyl Ethers (PBDE): < 1000 ppm
Bis(2-Ethylhexyl) phthalate (DEHP): < 1000 ppm	Benzyl butyl phthalate (BBP): < 1000 ppm
Dibutyl phthalate (DBP): < 1000 ppm	Diisobutyl phthalate (DIBP): < 1000 ppm

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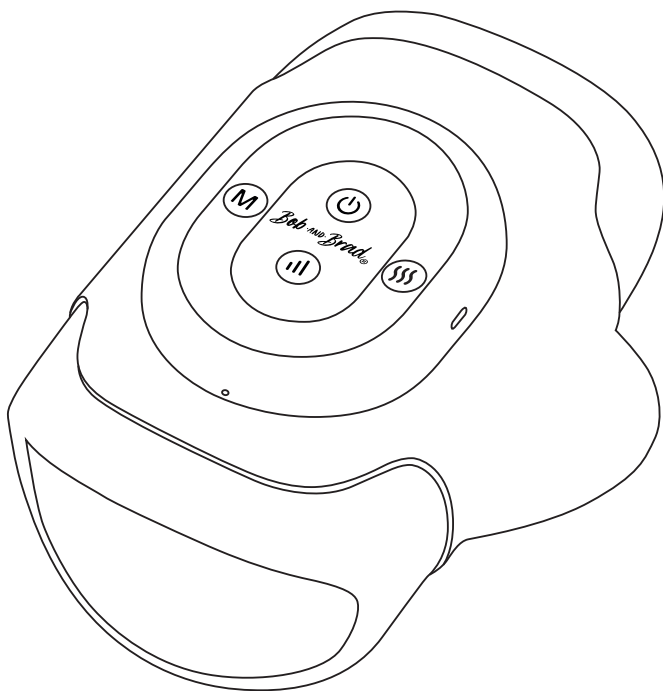Manufacturer:

Shenzhen Laige Technology Co.Ltd.

Room 401, 4/F, Building 5, Yunli Intelligent Park, Bantian Street,

Longgang District, Shenzhen, Guangdong, 518129, China

support@bobandbrad.com

EC Declaration of Conformity

We(UK Importer is responsible for this declaration

KOVA ASSOCIATES LTD

(Company name)

72a Kingston Road, Portsmouth, Hampshire, PO2 7PA, United Kingdom

(Company address)

KOVA.Compliance@outlook.com

(Email Address)

declare under our sole responsibility that the product

Tradename: BOB AND BRAD

Equipment: Hand Massager

Model No.: H50

ASIN: B0DTYCRSPN

to which this declaration relates is in conformity with the essential requirement set out in the Council Directive on the Approximation of the laws of the Member States relating to LVD Directive(2014/35/EU) & EMC Directive(2014/30/EU) & RoHS(2002/95/EC) product is responsible to affix CE marking, the following standards were applied:

EN 60335-1:2012+A11:2014+A13:2017

EN 60335-2-32:2003+A1:2008+A2:2015

EN 62233:2008

EN 55014-1:2017

EN 55014-2:2015

EN 61000-3-2:2014

EN 61000-3-3:2013

Full Name: Taozirui

Position: Manager

Signature: Taozirui

Date: 2021.5.6

EC Declaration of Conformity

We (EU Importer is responsible for this declaration)

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EN 61000-3-2:2014

EN61000-33:2013

Full Name: Taozirui Position: Manager

Signature: Taozirui Date: 2021.5.6

We Stand Behind Our Products!

2-Year Limited Warranty & Troubleshooting **Support@bobandbrad.com**

At Bob and Brad, your happiness comes first. Our products come with a **2-year** warranty, and our Customer Happiness Team is always here to help--anytime you need us. Happiness, Guaranteed.

Troubleshoot

Help you solve any issues with your device.

Product Knowledge

Answer questions about your device and best practices.

Initiate Returns

Provide instructions for replacements, returns and refunds.

Place New Orders

Help with placing orders online or via authorized resellers.