medisana®

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Heizkissen HP 515 / HP 516 / HP 517 / HP 518

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Device and controls

- ❶ Control unit HD20-A
- ❷ On/Off button
- ❸ Cable with plug
- ❹ Heating pad
- ❺ Plug connection
- ❻ Button to change heating level

Explanation of symbols

Do not use the device when it is folded!


Do not puncture the device!

Not suitable for children under 3 years!

The heating pad (without control unit) can be cleaned in a gentle wash cycle at max. 30°C!

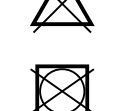
Do not bleach!

Do not dry in a tumble drier!

Do not iron!

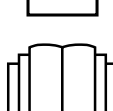Do not dry clean!

The device must only be used indoors!

IMPORTANT
Non-observance of these instructions can result in serious injury or damage to the device.

WARNING
These warning notes must be observed to prevent any injury to the user.

CAUTION
These notes must be observed to prevent any damage to the device.

NOTE
These notes give you useful additional information on the installation or operation.

Protection category II

LOT
LOT number

Manufacturer

Recycling symbols / codes:
These are used to provide information about the material and its proper use and recycling.

GB Important Instructions
Retain for future use

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- The device must only be used for its intended purpose as described in the instruction manual. Using the appliance for any other purpose invalidates the warranty.
- The device is not to be used for commercial use or medical purposes.
- When the device is likely to be used for an extended period, for instance by the user falling asleep, controls should be adjusted to a setting recommended for continuous use.
- Check the article carefully for signs of wear and/or damage before each use. Do not operate the device if you notice wear, damage or signs of improper use, or if it does not function correctly. Instead, contact an authorised service centre in such a case.

Safety Instructions – About the power supply

- Before connecting the heating pad to the power supply, please check that it is switched off and that the mains voltage of your power supply corresponds to the voltage listed on the label.
- Keep mains cable and control unit away from heat, humidity and liquids. Never touch the mains plug when your hands are wet or when you are standing in water.
- Only operate the heating pad with the corresponding switching unit HD20-A.
- Do not attempt to touch a heating pad which has fallen into water. Pull the plug out of the socket immediately.
- To disconnect the heating pad, always pull the mains plug out of the mains socket. Never pull on the cable!
- Do not carry, pull or twist the heating pad using the mains cable.
- Make sure that the power socket and plug are easily accessible so you can unplug the mains quickly if necessary.
- The cable and control unit of the appliance can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user must ensure (e.g. using cable ties) that all electrical cables are arranged in a safe way.

Safety Instructions – Specific groups of people

- The device complies with the guidelines on electromagnetic compatibility. If you wear a pacemaker, this could be affected. Consult your doctor and the manufacturer of your pacemaker before using this device.
- If you have any medical concerns, please discuss using the heating pad with your doctor before using it.
- Do not use the device on children or persons who are disabled, asleep or insensitive to heat (persons who cannot react on overheating).
- Children below the age of 3 years are not allowed to use this device, because they are not able to react on overheating.
- The device must not be used by young children over the age of 3 years, unless the control unit has been adjusted accordingly by its parents or legal guardians and unless the child has been well instructed how to use the control unit properly.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- The appliance can be used by children older than 3 years and younger than 8 years under supervision and with the control always set to minimum temperature value.
- Children shall not play with the appliance.
- Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!
- Do not use the heating pad on parts of the body which are swollen, inflamed or injured.
- Consult your doctor, if you have pain in muscles or joints for a prolonged period of time.

Safety Instructions – Operating the device

- Only use the heating pad as intended and in accordance with the instruction manual.
- The warranty will become void if it is used improperly.
- Do not use the unit for animals.
- This heating pad is intended only for private use and not for commercial purposes or for use in hospitals.
- The heating pad is subject to aging. Check heating pad, mains cable and control unit carefully for damage before each use. A faulty unit must not be used.
- Do not use the appliance if it is damaged, if its power cord is damaged, if it is not in perfect working order, or if it has been dropped or fallen into water. Send the unit to the service centre for repair to avoid potential hazards.
- Keep the unit out of rain and do not let it get wet. It must not be used out-of-doors.
- Never place the unit immediately next to a heater or radiator.
- Ensure that the unit does not come into contact with pointed or sharp objects and do not push pins into it.
- The heating pad must not be used if it has been folded, bent or crushed.
- Do not use the heating pad underneath quilts, blanket or pillows or cushions.
- Do not use the heating pad when wet and only use it in dry environments (not in the bathroom or similar environments).
- The control unit must not be placed on or under the heating pad during operation and it must not be covered.
- Do not sit on the heating pad, but place it like a cushion on the body part, as required.
- Do not leave the heating pad unattended. Always switch off the heating pad when not in use and pull the plug out of the socket. Never use a time switch.
- Do not fall asleep while the heating pad is switched on and the plug is inserted into the socket.
- Ensure that you follow the safety notes, as prolonged use at a high setting can result in burns to the skin. Check the skin areas under the heating pad regularly to prevent burns or blisters.
- If use of the heating pad is unpleasant or painful, discontinue using it immediately.

Safety Instructions – Maintenance and cleaning

- You may only clean the heating pad. To avoid potential hazards, you must never attempt to repair it yourself. Contact the service centre.
- Cleaning and user maintenance shall not be made by children without supervision.
- Before cleaning the unit, switch it off and unplug it from the mains. Never immerse the control unit in water or other fluids.
- Let the heating pad cool down completely before storing it. Store it in a place that is dry and free from dust, preferably in the original packaging, until you use it next.
- Do not crease the appliance by placing items on top of it during storage.
- If a fault occurs, do not attempt to repair the unit yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- A damaged power supply cable must be replaced with a power cable of the same type.

Scope of delivery

Please check first whether the unit is complete and does not show any damage.
If in doubt, do not operate the unit and contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Heating Pad with removable control unit HD20-A
- 1 Quick Start Guide (short instructions)

If you notice any transport damage when unpacking, please contact your dealer immediately.

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.

WARNING
Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Use

You can use the heating pad to gently heat and relax any part of the body. Use of the heating pad promotes blood circulation through sore muscles and increases wellbeing after a strenuous day. The heating pad reaches a pleasant temperature within a few minutes of switching on. The temperature can be adjusted as required.

Control Unit HD20-A

- ❷ On/Off button
- ❻ Button to change heating level
- ❸ LED display of the heating level (available levels: 1 ... 4)

Using the device

Make sure, the connection cable of the control unit is properly connected to the plug of the heating pad. Plug the unit into the mains outlet and press the On/Off button ❷. The device starts at lowest level (level 1), which is displayed at ❸.

After a few minutes you will be able to feel the pad getting noticeably warmer. If you require a higher temperature setting, press the Heat button ❻ several times until you reach the required level. (The maximum level is level 4, after that the level returns to the lowest level 1.) If the heating pad gets too hot for you, press the Heat button ❻ several times until you reach the required level. After a short time, the heating pad will cool down.

The heating pad may cool down after a short period of time. This is normal. The fast re-heating according to the setting will follow immediately.

For continuous use, you have to always select the lowest level 1. To switch the heating pad off, press the On/Off button ❷. The LEDs ❶ will go out, indicating that the unit is switched off. Unplug the unit from the mains outlet.

The heating pad will switch off automatically after approx. 90 minutes of continuous use (unless you have manually switched it off before). Unplug the unit from the mains outlet once you have finished using it.

Troubleshooting

If there is no function at all after switching on the heating pad (no heating level is displayed at ❸), please check if the mains plug is correctly plugged in the mains outlet. If this is the case, the heating pad is defect. Then please do not use the heating pad and contact the service center.

Cleaning and maintenance

- Before cleaning the pad, unplug it from the mains outlet and let it cool down for at least ten minutes.
- The heating pad is equipped with a removable control unit (HD20-A). Unplug the control unit from the heating pad by pulling out the connection cable from the plug on the heating pad.
- You may wash the heating pad gently by washing machine at 30°C.
- However, the best way is to put the heating pad in a bath with lukewarm water and some mild detergent and then gently squeeze it out. Rinse the heating pad several times to remove all the detergent.
- Never use strong detergents or hard brushes.
- To dry, lay the pad flat on a moisture-absorbent towel or mat.
- Only use the blanket again once it is completely dry.
- Reconnect the control unit to the heating pad in the correct manner.
- You should untangle the cable from time to time.
- Store the unit in a clean and dry place, laid out flat with nothing covering or resting on top of it.

Technical Data

Name and model:	medisana Heating Pad (with Control Unit HD20-A)
Power supply:	via Control Unit HD20-A input: 220–240V~ 50/60 Hz Length of the mains cable: approx. 1.8 m approx. 80 W (at max. heating level)
Heat output:	4 levels.
Heating levels:	The active level (1 ... 4) is displayed via four LEDs ❶.
Max. operating time:	approx. 90 minutes, with automatic switch-off
Dimensions (L×W×H):	approx. 40 × 30 × 0.5 cm
Weight:	approx. 500 g
Operating conditions:	Use only in dry indoor areas
Storage conditions:	Laid out flat and dry

Name	Colour	Item number / EAN number:
HP 515	denim blue	61130 / 40 15588 61130 8
HP 516	light grey	61131 / 40 15588 61131 5
HP 517	savage green	61132 / 40 15588 61132 2
HP 518	light purple	61133 / 40 15588 61133 9

Disposal

 This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

- The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:
 - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - Accessories which are subject to normal wear and tear.
- Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH, Carl-Schurz-Str. 2, 41460 NEUSS, GER-MANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under <https://docs.medisana.com/61130>.
Information about service can be found here: www.medisana.com/servicepartners

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medisana GmbH , Carl-Schurz-Str. 2, 41460 NEUSS, ALMANYA