


Product brochure

08.08.2026

KERN ADT 100-4

Analytical balance with flexible functionality and touch display

ENGLISH

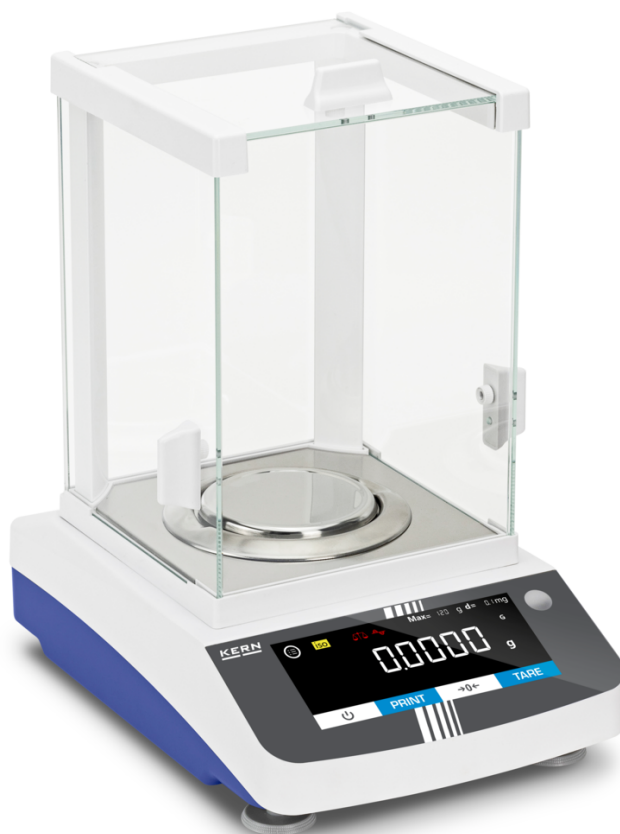

CONTACT

KERN & SOHN GmbH
Ziegelei 1
72336 Balingen
Germany

+49 7433 9933-0
info@kern-sohn.com
www.kern-sohn.com



APPLICATION EXAMPLE



FEATURES

- The wide range of functions of this range of analytical balances means that they qualify for use in the pharmaceutical industry as well as other fields. These include, for example, the statistical function, weighing with tolerance range, counting function, percentage weighing, totalizing function and many more
- Modern touch display with convenient operation enables you to, for example, switch the units or start adjustment straight from the display
- The Windows Direct function allows weights to be transferred directly from the balance to a Windows application via a USB device connection, without the need for manual entry
- Draught shield standard
- Specialised functions, such as the statistical function make it easier to use efficiently in industry and the laboratory
- GLP/ISO record keeping with date, time and identification number of balance adjustment or weighing process
- Assists in data integrity in accordance with U.S. FDA 21 Part 11 (for example weighing result, sample ID, user name, scales ID, ...)
- Data query and remote control of the balance using a computer or CRM/ERP systems using the KERN Communication Protocol
- RS-232 data interface and USB (Device) for the transfer of weighing data
- Menu language EN

TECHNICAL DATA

- Backlit LCD display, digit height 20 mm
- Dimensions weighing surface, stainless steel, Ø 90 mm
- Overall dimensions W×D×H 207×318×360 mm
- Net weight approx. 6 kg
- Permissible ambient temperature 15 °C/25 °C

MAIN SCOPE OF APPLICATIONS

- Laboratories
- Microanalysis
- Pharmaceutical applications
- Material testing
- Precision engineering
- Research and development
- Chemical industry
- Food industry
- Plastics industry
- Quality laboratories

PICTOGRAMS

Standard

Option



MODEL OVERVIEW

Model	Weighing capacity [Max]	Readability [d]	Reproducibility	Linearity	Dimensions weighing surface
ADT 100-4	120 g	0,0001 g	0,0003 g	± 0,0003 g	90 mm
ADT 200-4	220 g	0,0001 g	0,0003 g	± 0,0003 g	90 mm
ADT 300-4	320 g	0,0001 g	0,0004 g	± 0,0004 g	90 mm

ACCESSORIES

Model	Description
315-107-103	Safety set 315-107-103, included weight 5 g (E2); 100 g (E2)
316-07	Test weight 316-07 (OIML Class E2, Nominal value 100 g, Compact form, stainless steel polished)

Model	Description
572-926	RS232 interface cable 572-926
ABP-A01	Ioniser ABP-A01
ABS-A08	Protective dust cover ABS-A08
YBA-A24S05	Screen protector YBA-A24S05
YBI-01A	Ioniser YBI-01A
YDB-03	Density determination set YDB-03
YKC-01	Thermal printer YKC-01
YKG-01	Dot Matrix Printer YKG-01
YKH-01	Thermal printer YKH-01
YKI-10	RS232-Ethernet adapter YKI-10 (RS-232 standard, Ethernet, RS-485)
YKI-11	RS232-Ethernet adapter YKI-11 (Ethernet, RS-232 standard - 2x)
YKI-12	RS232-WiFi adapter YKI-12
YKT-01	Needle printer YKT-01
YPS-03	Antivibration table YPS-03
YPS-04	Anti-vibration plate KERN YPS-04 (400×450×60 mm, 30 kg)
YPS-05	Anti-vibration plate KERN YPS-05 (565×450×60 mm, 42 kg)

SOFTWARE

Model	Description
SCD-4.0-DL	Software SCD-4.0-DL (Scope of supplies: Download link for 1 license)
SCD-4.0-DLS05	Software SCD-4.0-DLS05 (Scope of supplies: Download link for 5 licenses)
SCD-4.0-PRO-DL	Software SCD-4.0-PRO-DL (Scope of supplies: Download link for 1 license)
SCD-4.0-PRO-DL-AS05	Software SCD-4.0-PRO-DL-AS05 (Scope of supplies: Download link for 5 licenses)

SERVICES

Model	Description
963-101	DAkkS calibration certificate 963-101
965-411	Assembly of components 965-411 (< 50kg)
969-517	Confirmation of the manufacturer's specification (only in combination with DAkkS calibration certificate) 969-517
970-108	Warranty extension (+2 years) 970-108

KERN Pictogram Overview

	Internal adjustment: Adjusting accuracy using internal adjustment weight
	Data interface RS-232: For connection to printers, PCs or other peripheral devices for data transfer over longer distances. Network in bus topology possible
	Data Interface USB: For connection to printers, PCs or other peripheral devices
	KERN Communication Protocol (KCP): a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems
	GLP/ISO protocol: The measuring device generates a GLP-compliant printout
	Value & time: The measuring device displays the value, date, and time, regardless of the connected printer
	Piece counting: Reference quantities selectable. Display can be switched from piece to weight
	Net total: The weight values of the recipe ingredients can be added up and the total weight of the recipe printed out
	Density determination: The density of liquids and solids with a density of $\leq / \geq 1$ is determined directly in the measuring device
	Percentage determination: determining the deviation in % from the target value (100 %)
	Units can be switched to e.g. non-metric units
	Measuring with tolerance range (checkweighing): Upper and lower limit values can be programmed, e.g. for sorting, portioning or for measuring coatings or material thicknesses. The process is supported by an acoustic and/or optical signal
	Hold function: When the weighing conditions are unstable, a stable weight is calculated as an average value
	Suspended weighing: Load support with hook on the underside of the balance
	Universal plug-in power supply: With universal input and optional input socket adapters for EU, CH, GB or EU, CH, GB, US or EU, CH, GB, US, AUS
	Weighing principle electromagnetic force compensation: Coil inside a permanent magnet. For the most accurate weighings
	Package shipment by courier service: The duration of internal product provision is one day
	Accredited calibration (DKD): The duration of the accredited calibration is 3 working days