VEVOR

Affordable. Reliable. Home Improvement.

MICROSCOPES

MODEL: H-20A, H-20B

https://youtu.be/o9T8qo1gGlQ

MODEL: H-20A, H-20B

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

⚠ WARNING:

Please read this manual carefully before using the product. Failure to do so may result in serious injury. SAVE THESE INSTRUCTION

1. Never, under any circumstances, look directly at the Sun, another bright source of light or at a laser through this device, as this may cause permanent retinal damage and may lead to blindness
2. Do not use the microscope in extremely cold or hot and humid conditions. Store the device in a dry, cool place away from hazardous acids and other chemicals, away from heaters, open fire and other sources of high temperatures.
3. This microscope is not water proof, please avoid splashing water or any other liquid directly. In the event any liquid or solution has leaked into the microscope first unplug the power cord and refer to professional service personnel. Does not contain any self serviceable parts.
4. Do not use the device for lengthy periods of time, or leave it unattended in direct sunlight. When using the microscope, try not to use it near flammable materials or substances (benzene paper, cardboard, plastic, etc.), as the base may heat up during use, and might become a fire hazard.
5. Do not touch the optical surfaces with your fingers. Do not use any corrosive or acetone based fluids to clean the optics. Do not attempt to dismantle the microscope and refer to professional service personnel.
6. Protect the device from sudden impact and excessive mechanical force. Do not apply excessive pressure when adjusting focus. Do not over tighten the locking screws.
7. Battery and main power may not be used at the same time. Unplug the power cord before using battery. Always use the power supply with the proper voltage.
8. Grounding the microscope is recommended. Electric shock may damage sensitive semi conductors inside the microscope. And avoid handling the power cord with wet hands.
9. Operate the microscope as directed in the manual and extreme caution should be practiced at all times. Children should use the device under adult

supervision only. Seek medical advice immediately if a small part or a battery is swallowed.

MAINTENANCE

1. The microscope should be kept out of direct sunlight and stored in a cool, dry place, free from dust, fumes and moisture. It should be stored in a case or any place covered with a hood to protect it from dust.
2. The microscope has been carefully tested and inspected. Since all lenses have been carefully aligned, they should not be disassembled. If any dust has settled on the lenses, blow it off with an air blower or wipe off with a clean soft came hairbrush. In cleaning mechanical parts and applying non-corrosive lubricant, take special care not to touch the optical elements, especially the objective lenses.
3. If you are not using your microscope for extended periods of time, store the objective lenses and eyepieces separately from the microscope. When disassembling the microscope for storage, always put the covers on the Exchangeable Objective Lenses opening to prevent dust settling inside the lenses. Also keep the neck of the head covered.
4. Always unplug the microscope from a power source before opening the base or changing the illumination lamp. Regardless of the lamp type (halogen or incandescent), give it sometime to cool down before trying to change it, and always change it to a lamp of the same type.


TECHNICAL SPECIFICATIONS

Model	H-20A, H-20B
Adapter Input	AC110~240V,50/60Hz
Rating(s)	DC 5V, 2A
Product Dimensions (LxWxH)(mm)	135*170*300
Net Weight(kg)	About 2.04

*Products such as specifications, appearance, and design are subject to modification without prior notice.

1) Eyepices: WF10X/WF25X, 2X magnifier
2) Objective Lens: 4x,10x,40x
3) The total magnification of the microscope = (the magnification of the objective lens in use) X (the magnification of the eyepiece in use) X (the magnification of the extender)
4) Light Source: LED light
5) Power Supply: AC Adapter

COMPONENT NAME

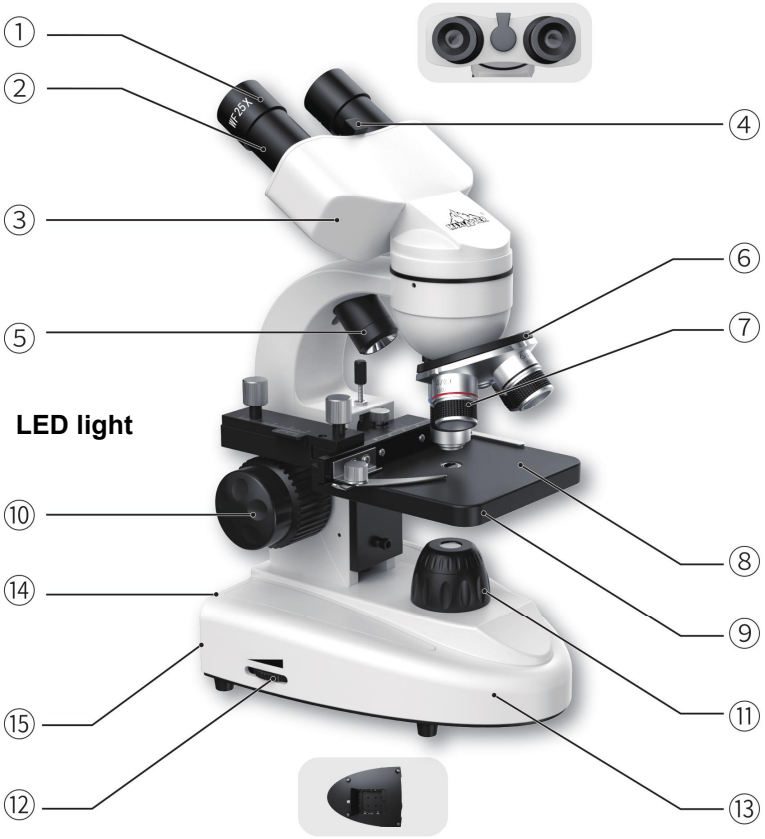① Eyepiece
② Eyepiece Tube
③ Binocular Head
④ Interpupillary Distance Adjuster
⑤ Upper LED Light
⑥ Nose-piece
⑦ Objective Lens
⑧ Stage
⑨ Diaphragm
⑩ Coarse/Fine Focus Knob
⑪ Bottom LED Light
⑫ Brightness Control Knob
⑬ Base
⑭ Battery Bin
⑮ Power Interface

Specifcations:

1. Observing Binocular Head: 45° inclined eyepiece tube 360° rotation head
2. Mechanical Stage: 90*90mm stage with scale and clip
3. Focus Knob: Coarse and Fine focusing knob

INSTRUCTIONS FOR USE OF MICROSCOPE

STEP1
Connect the power supply


STEP2
Insert the eyepiece into the eyepiece tube. Use the interpupillary distance adjuster to adjust the distance between the eyepieces.

STEP3
Place the microscope slide on the stage and secure it with a stage clip.

STEP4

Turn the nosepiece to bring the low-magnification objective lens (e.g. 4x) to the working position.

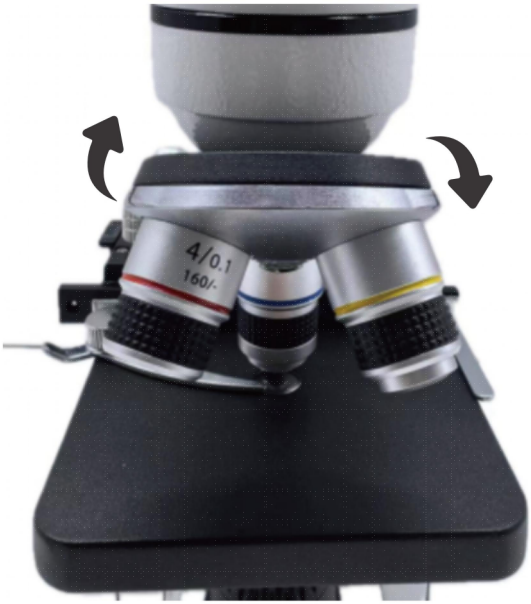
STEP5

Turn on the upper and bottom LED lighting.

STEP6

Turn the moving ruler to move the specimen slide so that the center of the slide is close to the center of the light hole.

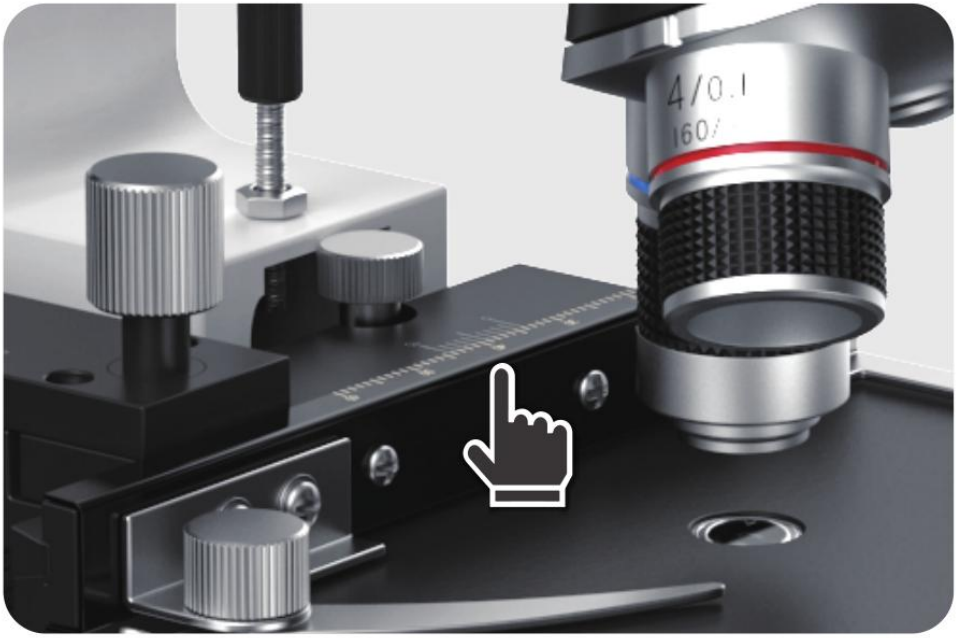STEP7

Use the coarse focus screw to adjust the focal length of the objective lens until the image of the specimen is bright and clear.

STEP8

Use the fine focus screw to finely focus the image to obtain a satisfactory clear image.

STEP9

Turn the nosepiece to change other objective lenses for observation, and then readjust the thickness of the focus screw to obtain the clearest image

TROUBLESHOOTING GUIDE

Instrument Part

Trouble	Cause	Solution
Dark edges or uneven light and shade in the field of view	The nosepiece is not in the position	The nosepiece is turned to the positioning position
There is dirt in the field of view	There is dirt on the lens (objective, eyepiece, condenser lens) or specimen slide	Wipe clean
The image quality is very poor (low resolution, poor contrast)	There is no cover slip on the specimen	Cover slip on specimen
	Coverslip is too thick or too thin	Coverslips with standard 0.17mm thickness
	Immersion oil on dry objectives (especially 40X)	Wipe clean

Operation Part

Trouble	Cause	Solution
The objective lens can easily break into pieces during operation.	Platform limit is too high	Lower the platform limit screw
The image cannot be focused when observed with a high-magnification objective lens	The specimen piece is placed upside down or the coverslip is too thick.	Turn the specimen slide over and cover slip with a standard thickness (0.17mm)
Contact the glass when the objective lens is switched from low	Coverslip is too thick	Use standard thickness coverslips (0.17mm)

magnification to high magnification		
The specimen does not move smoothly	The specimen holder is not tightened	Retighten

Electrical Part

Trouble	Cause	Solution
The bulb does not light up when the switch is turned on	No power	Check the connection of wires and the presence of power
	The light bulb is broken	Replace light bulb
The lighting is not bright enough	Power supply voltage too low	Increase voltage
The light bulb flickers or the brightness is unstable	The light bulb is about to break	Replace

Correct Disposal:

This product is subject to the provision of european Directive 2012/19/EU. The symbol showing a wheelie bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to acollection point for recycling electrical and electronic devices.

FCC Information:

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment!

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

1) This product may cause harmful interference.

2) This product must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this product not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

Note: This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This product generates, uses and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception,

which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

· Reorient or relocate the receiving antenna.

· Increase the distance between the product and receiver.

· Connect the product to an outlet on a circuit different from that to which the receiver is connected.

· Consult the dealer or an experienced radio/TV technician for assistance.

IMPORTANT!

Avoid touching the top of the machine while it's in use, as it can become very hot. There commended intermittent operation cycle for this machine is 30 minutes on and 10minutes off. However, we advise allowing the machine to cool down for at least 15minutes after 30 minutes of use. Never place items on top of the machine under any circumstances.

Manufacturer: Shanghaimuxinmuyeyouxiangongsi
Address: Shuangchenglu 803nong11hao1602A-1609shi, baoshanqu, shanghai 200000 CN.
Imported to AUS: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australia
Imported to USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730

UK	REP

YH CONSULTING LIMITED. C/O YH Consulting Limited Office 147, Centurion House, London Road, Staines-upon-Thames, Surrey, TW18 4AX

EC	REP

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Importerad till Australien: SIHAO PTY LTD. 1 ROKEVA STREETEASTWOOD NSW 2122 Australien
Importerad till USA: Sanven Technology Ltd. Suite 250, 9166 Anaheim Place, Rancho Cucamonga, CA 91730