
medisana®

51179 BU 535 Voice 6gs 31-Jul-2025 Ver. 1.7

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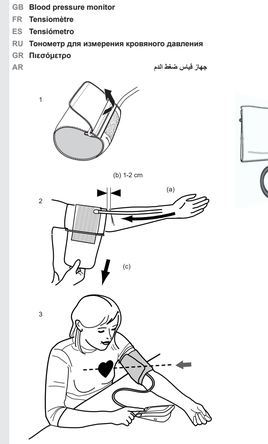medisana GmbH,
Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

ANDON HEALTH Co., Ltd.
No. 3 Jingjing Street, YaAN Road, Nankai District, Tianjin 300190, CHINA

EC REP **iHealthLabs Europe SAS**
36 rue du Ponthieu, 75008, Paris, France

Name: BU 535 Voice (Model: KD-5920V)
REF **51179**
EAN **4015588 51179 0**

Manual Service

MD CE 0197

GB

IMPORTANT INFORMATION!
KEEP IN A SAFE PLACE!

Please read this instruction manual carefully, in particular the safety instructions, before using the device. Keep the instruction manual in a safe place for later reference. If you pass the device on to a third party, this instruction manual must remain with the device.

- The unit is intended only for use in the home. Consult your doctor prior to using the blood pressure monitor if you are concerned about health matters.
- Only use the unit for its intended purpose in accordance with the instructions for use. Warranty claims become void if the unit is misused.
- Please consult your doctor prior to using the unit if you suffer from illnesses such as arterial occlusive disease.
- The unit may not be used to check the heart rate of a pacemaker.
- Like on any other oscillometric blood pressure measurement devices, certain medical conditions can affect the measurement accuracy, among others: cardiac arrhythmias, low blood pressure, circulatory disorders, states of shock, diabetes, pregnancy, pre-eclampsy etc. Always consult a physician before you use the device.
- This device is not suitable to provide a medical diagnosis. The measurements are for informational purposes only to get information about treatment and medicine, which are necessary for you.
- The detection of an irregular heartbeat is not a substitute for a professional heart examination. Check with your doctor if any treatment or medical action is necessary for you.
- Irregularity of pulse or arrhythmia can lead to difficulties in recording a correct reading when measurements are taken using oscillometric blood pressure devices. This device is electronically equipped to detect irregular heartbeat which it indicates with the symbol ♡. In this case, contact your physician.
- Please contact your physician about the specific risks of cuff pressure in your individual case.

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SN LOT MD EC REP IP UDI

SAFETY NOTES FOR BATTERIES
- Do not disassemble batteries.
- Replace all batteries if the unit display shows the low battery symbol.
- Never leave any low battery in the battery compartment since it may leak and cause damage to the unit.
- Increased risk of leakage! Avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with any of these parts, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed, seek medical attention immediately!
- Replace all of the batteries simultaneously!
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones!
- Insert the batteries correctly; observing the polarity!
- Remove the batteries from the unit if it is going to be used for an extended period of 3 months or more!
- Keep batteries out of children's reach!
- Do not short circuit batteries! Do not short circuit! Do not throw into a fire! **There is a danger of explosion!**
- Do not throw used batteries into the household refuse; put them in a battery collection point, at the shop where they were purchased!

Proper use
This fully automatic electronic blood pressure monitor is intended for measuring the blood pressure at home. It is a non-invasive blood pressure measurement system to measure the diastolic and systolic blood pressure and pulse of an adult using a non-invasive technique by means of a cuff, which needs to be fitted on the upper arm. The cuff size is suitable for arm circumferences between 22 - 42 cm.
Intended user: Medical professionals or lay person.
Contraindications: The device is not suitable for measuring the blood pressure of children. Ask your doctor before using it on older children.
This device is not suitable for persons who suffer from strong arrhythmia.

Unit and LC-Display

❶ Upper arm cuff	❼ Pump	⓫ LC-display	⓯ ▲ / ▼ button (START/STOP)
❷ M - button (Memory Recall)	❽ Battery compartment (on underside)	⓬ Push-in Connector for air tube	⓰ Display of time
❸ Air tube	❾ Display of Systolic Value	⓭ Display of Diastolic Pressure	⓱ Display of Pulse Rate
❹ Pulse symbol	❿ Symbol for zero pressure search	⓮ Memory Location number	⓲ M - symbol (memory)
❺ Blood Pressure Indicator	⓫ Change Battery symbol		⓳ Unit has detected irregular heartbeat

Scope of delivery
The delivery includes:
- 1 **medisana** Blood pressure monitor **BU 535 Voice**
- 1 storage pouch
- 4 batteries (type AAA, LR03) 1.5V
- 1 instruction manual
- 1 cuff with air tube

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your supplier or your service centre.
If you notice any transport damage during unpacking, please contact your supplier without delay.

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required.

WARNING
Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

How is blood pressure measured?
The **medisana BU 535 Voice** is used to measure blood pressure at the upper arm. The measurement is carried out by a microprocessor, which, via a pressure sensor, measures the vibrations resulting from the inflation and deflation of the cuff over the artery.

Blood pressure classification

systolic mmHg	diastolic mmHg	Blood pressure indicator ❷ / ❼
≥ 180	≥ 110	severe hypertension — red
160 - 179	100 - 109	moderate hypertension — orange
140 - 159	90 - 99	mild hypertension — yellow
130 - 139	85 - 89	high-normal blood pressure — green
120 - 129	80 - 84	normal blood pressure — green
< 120	< 80	optimal blood pressure — green

WARNING
Blood pressure that is too low represents just as great a health risk as blood pressure that is too high! Fits of dizziness may lead to dangerous situations arising (e.g. on stairs or in traffic)!

Influencing and evaluating readings
- Measure your blood pressure several times, then record and compare the results. Do not draw any conclusions from a single reading.
- Your blood pressure readings should always be evaluated by a doctor who is also familiar with your personal medical history. When using the unit regularly and recording the values for your doctor, you should visit the doctor from time to time to keep him updated.
- If the systolic and diastolic readings seem unusual (too high or too low) on several consecutive occasions, despite correct use of the unit, please inform your doctor. This also applies to the rare occasions when an irregular or very weak pulse prevents you from taking readings.

Starting up: Insert / replace batteries
You must insert the batteries provided before you can use your unit. The lid of the battery compartment ❽ is located on the underside of the unit. Open it by pushing the cover outward and remove it. Insert the 4 AAA LR03 type 1.5 V batteries supplied. Ensure correct polarity when inserting the batteries (as marked on the battery compartment cover). Close the battery compartment. Replace the batteries immediately, if the change battery symbol 🔋 appears on the display 3 or if nothing appears on the display after the unit has been switched on.

Settings: Display format, date, time
For correct operation of the device it is necessary to set time and date. Immediately after inserting the batteries the device is in setting mode for 1 minute. Alternatively, it will switch to standby-mode. You may enter the setting mode at any time by pressing and holding the ◔ - button ❻ and the M - button ❷ at the same time for approx. 2 seconds. In setting mode, firstly the clock format (12 or 24 hours format) can be set. The value can be adjusted with the M - button ❷. Press the ◔ - button ❻ afterwards to enter the date and time setting mode. Adjust the year, month, day, hour and minute figure in sequence. You may adjust the flashing value by pressing M - button ❷ and confirm the value and go to next step by pressing the ◔ - button ❻. After the minutes have been set and ◔ - button ❻ is pressed, all adjustments have been done and the device will switch itself to standby mode again. If no button is pressed during approx. 1 minute during the setting process, the device will automatically switch to standby-mode. If the batteries are changed, the set information is lost and needs to be entered again.

Setting the voice language
The **medisana** Blood pressure monitor **BU 535 Voice** features a voice output.
To adjust a language, press and hold the ▲ button ⓯ while the device is switched off. „L01", „L1" or similar appears in the display - the „0" represents a deactivated voice output and the other figures represent a language as follows:

Art. 51179:

1 = English	2 = French	3 = Arabic
4 = Russian	5 = Spanish	6 = Greek

Hold the **M** - button ❷ pressed down, until the desired language is set and then release the button. After the language is selected, you may use the **M** - button ❷ again in the same way to adjust the output volume.

Fitting the cuff
- Push the short piece of the air tube into the hole ❸ on the left side of the unit prior to use.
- Slide the open end of the cuff through the metal bracket so that the Velcro fastener is on the outside and it becomes a cylindrical form (Fig.1). Slide the cuff over your left upper arm.
- Position the air hose in the middle of your arm in line with your middle finger (Fig.2) (a). The lower edge of the cuff should be 2 - 3 cm above the crease of the elbow. (b). Pull the cuff tight and close the Velcro fastener (c).
- Measure the pulse on your bare arm.
- Only position the cuff on the right arm if it cannot be used on the left arm. Always carry out measurements on the same arm.
- Correct measuring position for sitting (Fig.3).

Taking a blood pressure measurement
After the cuff has been appropriately positioned the measurement can begin.
- Switch the unit on by pressing the ◔ - button ❻. The display will shortly show all available symbols.
- The **medisana** blood pressure monitor **BU 535 Voice** gives you the opportunity to assign readings to one of two user memories, each of which provides 120 memory slots. After switching on the device, press the M - button ❷, to choose between user 1 and user 2. If no selection is made within 5 seconds, the following measurement will be assigned to the current shown user.
- The last value measured for this user is shown in the display (if no value is available „0" is shown). The device now starts to deflate the cuff completely - the symbol for zero pressure search ⓾ is shown in the display.
- The unit is now ready for measurement and it automatically inflates the cuff slowly in order to measure your blood pressure. The rising pressure in the cuff is shown on the display. The unit inflates the cuff until sufficient pressure has built up for a measurement. Then the unit slowly releases air from the cuff and carries out the measurement. When the device detects the signal, the heart symbol ♥ ❹ starts to flash on the display.

WARNING
Do not take any therapeutic measures on the basis of a self measurement. Never alter the dose of a medicine prescribed by a doctor.

- When the measurement is finished, the cuff is deflated completely. The systolic and diastolic blood pressure and the pulse value appear on the display ❶. The blood pressure indicator ❷ is shown next to the relevant coloured bar depending on the blood pressure classification. If the unit has detected an irregular heartbeat, the symbol ♡ ⓳ will also be shown.
- The readings are automatically saved in the selected memory (1 / 2). Up to 120 measured values with date and time can be stored in each memory along with time and date.
- The measurement readings remain on the display. The unit switches off automatically after approx. 1 minute if no buttons are pressed. The unit may also be switched off using the ◔ - button ❻.

Discontinuing a measurement
If it is necessary to interrupt a blood pressure measurement for any reason (e.g. the patient feels unwell) the ◔ - button ❻ can be pressed at any time. The device immediately decrease the cuff pressure automatically.

Display stored values
This unit features 2 separate memories, each with a capacity of 120 memory slots. Results are automatically stored in the memory selected by the user. To call up the measured values stored, press the M - button ❷. Press the M - button ❷ when the device is switched off to activate the memory recall mode. A user memory (1 / 2) along with the number of the memory location is shown. Change the user memory with the ◔ - button ❻. Press the M - button ❷ to show the average values of the last 3 measurements of this user (along with „A3"). With each press on the M - button ❷ the display proceeds to show the average values for all measurements of this user memory ("AL"), the average value of the measurement results in the last 7 days taken between 5 o'clock and 9 o'clock in the morning ("AM") and the average value of the measurement results of the last 7 days taken between 6 o'clock and 8 o'clock in the evening ("PM"). Now press the M - button ❷ again to display the last measurement value of the user. Repeatedly pressing the M - button ❷ displays the respective values measured previously. If no further values are saved, the display will show a „0". If no button is pressed, the unit automatically switches itself off after approximately 1 minute. Memory recall mode may also be exited at any time by pressing the ◔ - button ❻, which will also switch the unit off. The oldest value is deleted if 120 measured values have been stored in the memory and a new value is saved.

Delete saved values
If you want to delete the saved measurements of a user, first switch the unit into memory recall mode (see „Display stored values" to display any record of the user. Now press and hold the M - button ❷ for approx. 3 seconds, The display will show „CL" and „00" and the values are deleted.

Error messages and error remedying
The following clarifications in the display indicate problems or errors:

Display	Cause and correction
Hi or Lo	The measurement result is out of the measurable range. Repeat the measurement and contact your physician if necessary.
Unusual measurement results	Cuff not fitted properly. Movement or speaking during measurement; Incorrect body posture; Heart arrhythmia. Repeat the measurement and contact your physician if necessary.
	Replace all batteries.
Er 0 / Er 1 / Er 2	Pressure system instable at start of the measurement / Error in measuring the systolic pressure / Error in measuring the diastolic pressure. **Do not move!**
Er 3 / Er 4	Pressure system is blocked or sleeve is too tight during inflation / Pressure system has a leak or sleeve is too loose during inflation. **Put the sleeve on correctly!**
Er 5 / Er 6 / Er 7 / Er 8 / Er A	Sleeve pressure above 300 mmHg / Longer than 160 seconds with sleeve pressure over 15 mmHg / EEPROM access error / Error when testing the functions on the device / Failure of pressure sensor. **Perform measurement again if errors occur.**
No reaction after pressing button or insertion of batteries	Remove the batteries, wait for 5 minutes and reinsert the batteries.

Please contact the manufacturer resp. the service centre if you can't solve the problem. Do not disassemble the unit by yourself.

Cleaning and care
- Remove the batteries before cleaning the unit. Clean the unit and the cuff using a soft cloth lightly moistened with a mild soapy solution. Never use abrasive cleaning agents, alcohol, naphta, thinner or gasoline etc. Never immerse the unit or any component in water. Be cautious not to get any moisture in the main unit. Do not use the unit again until it is completely dry. Only inflate the cuff when it is around your arm. Be careful not to bend, squeeze or damage the air tube. Store the unit in a clean and dry location.

Servicing and calibration: The unit has been calibrated by the manufacturer for a period of two years' use. In the case of commercial use, the unit must be serviced and recalibrated at least once every two years. This calibration will be charged for and can be carried out by an appropriate authority or an authorised service centre - in accordance with the regulations for the operation of medical equipment. The device is designed for maintaining the safety and performance characteristics for a minimum of 10,000 measurements resp. 5 years. The cuff integrity is maintained after 1,000 open/close cycles of the closure.

Hints: 1) The device shall not be serviced by the user.
2) We can provide product circuit diagrams and repairable component information to qualified maintenance service personnel if necessary.

Guidelines and standards
The device complies with all fulfils the requirements of the EU standard for non-invasive blood pressure monitors. It is certified in accordance with the EC mark (conformity mark) 'CE 0197'. The blood pressure monitor complies with the European standards EN 1060-3: 1997 + A2: 2009 ; ISO81060-2:2018.

Technical specifications

Type designation	: **medisana** Blood pressure monitor **BU 535 Voice**, not suitable for use in an oxygen-rich environment
Display system	: Digital display
Memory slots	: 2 x 120 measurement values
Measuring method	: Oscillometric
Power supply	: ⎓ 4 x 1.5 V batteries AAA LR03
Blood pressure measuring range	: 40 – 260 mmHg; Systolic: 60-260mmHg, Diastolic: 40-199mmHg, Cuff pressure: 0-300mmHg
Pulse measuring range	: 40 – 180 beats/min.
Maximum error tolerance for static pressure	: ± 3 mmHg
Maximum error tolerance for pulse rate	: Less than 60 beats/min: ±3 beats/min; More than 60 (incl.) beats/min: ±5%
Pressure generation	: Automatic with pump
Deflation	: Automatic
Automatic switch-off	: After approx. 1 minute
Operating conditions	: +10 °C to +40 °C; ≤ 85 % relative humidity; Air pressure 80 kPa - 105 kPa
Storage and transport conditions	: -20 °C to +50 °C; ≤ 85 % relative humidity; Air pressure 80 kPa - 105 kPa
Dimensions (main unit)	: approx. 150 x 95 x 41 mm (L x W x H)
Cuff	: 22 - 36 cm for adults
Weight	: approx. 235g without batteries
Battery life	: approx. 700 times
Replacement parts	
	- Cuff, L: 22 - 30 cm for adults with average upper arm circumference, Art. No. 51135 / EAN 40 15588 51135 6
	- Cuff, L: 30 - 42 cm for adults with upper arm of large circumference, Art. No. 51126 / EAN 40 15588 51126 4

Disposal — Waste Electrical and Electronic Equipment (WEEE)
This product is subject to European Directive 2012/19/EU on waste electrical and electronic equipment and is marked accordingly. Never dispose of electronic devices with household waste. Please seek information on local collection points and regulations with regard to the correct disposal of electrical devices. Correct disposal helps to protect the environment and human health.
All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your local authority or your supplier for information about disposal.

Warranty and repair terms
Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.
The following warranty terms apply:
1. The warranty period for **medisana** products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 c. Damage to the unit has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 d. Applied parts which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

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